

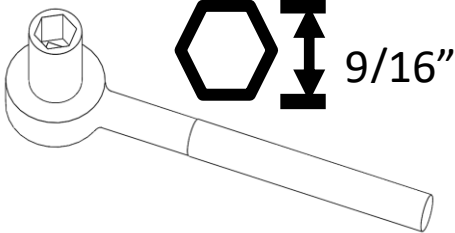
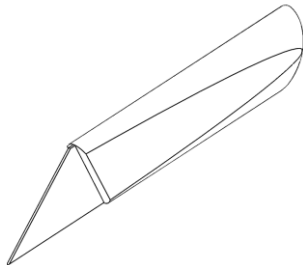
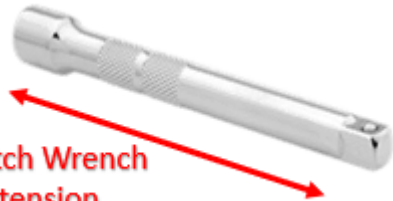



ENGLISH AND METRIC ALLEN WRENCHES

 9/16"	
 3/4 OPEN END WRENCH	 6" Ratch Wrench Extension
	Report any shipping Damage to IDI at 978-829-0009

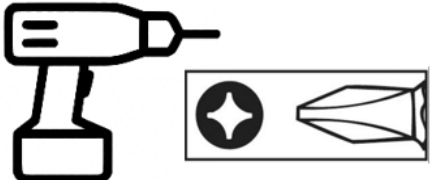


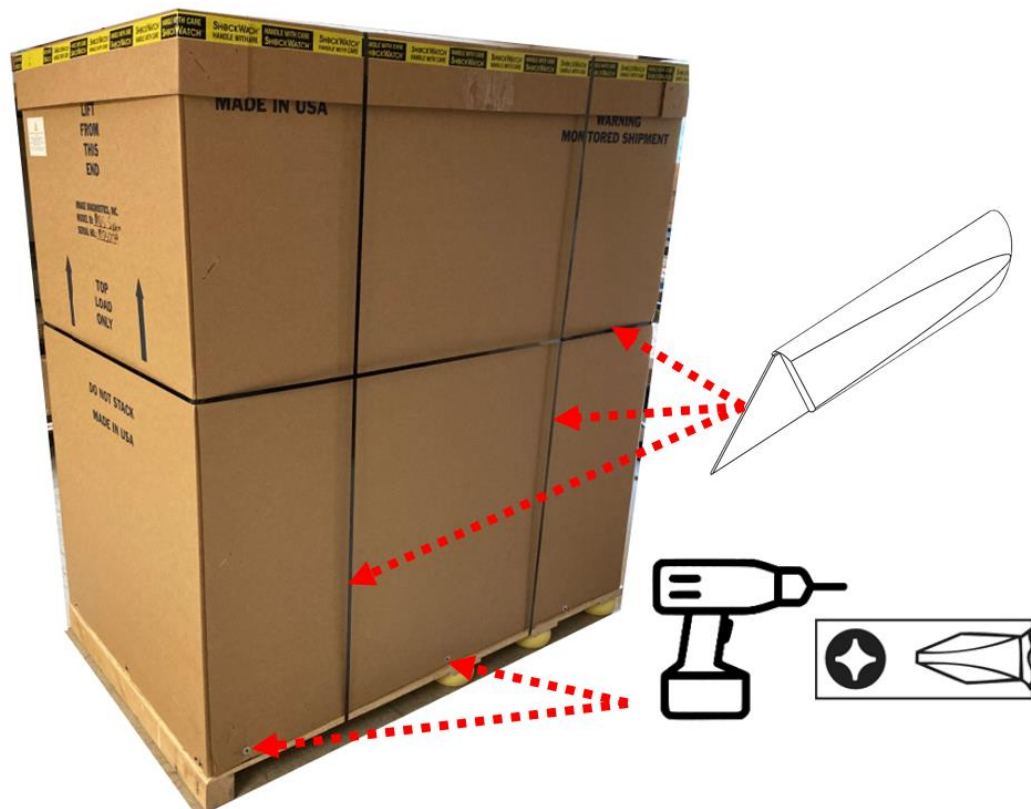

1.

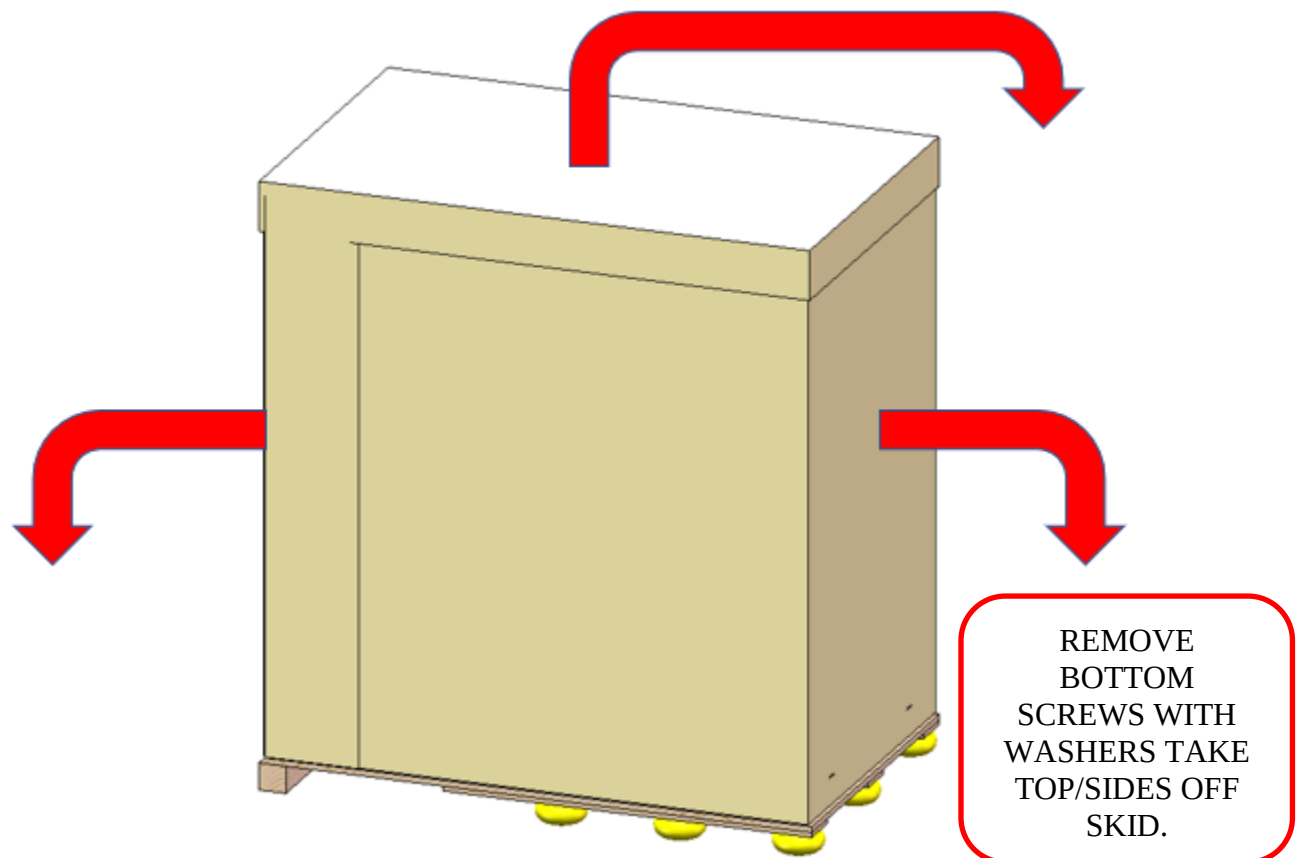


2.

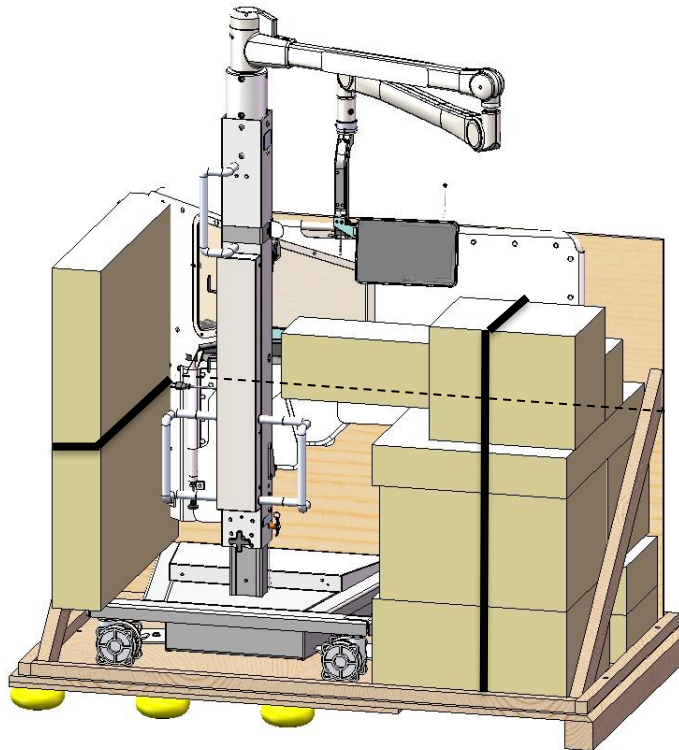


CUT ALL STRAP
AREAS,
MOVE BOXES
OFF SKID, MOVE
BOXED
COMPONENTS
TO DESIGNATED
LOCATION FOR
IDI REP. TO
UNBOX

3.

BARRIER: REMOVE
ALLENS, ROTATE WHITE
BLOCK OUT, MOVE
BARRIER PIPE INTO
WHITE BLOCK SLOT,
COLUMN STRAP AROUND
ARM EXTENSION

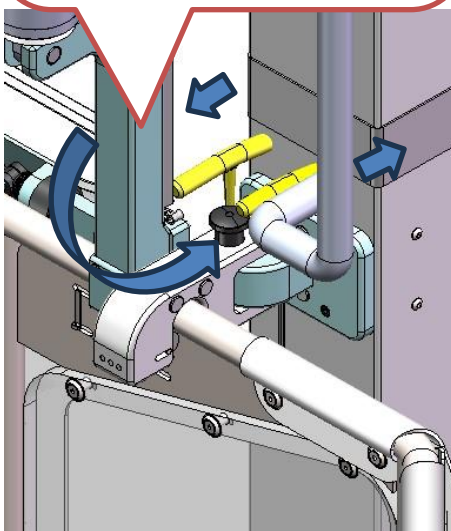


REMOVE
ANGLED 2X4
STUDS, USING
TOOLS, REMOVE
STEEL PLATES

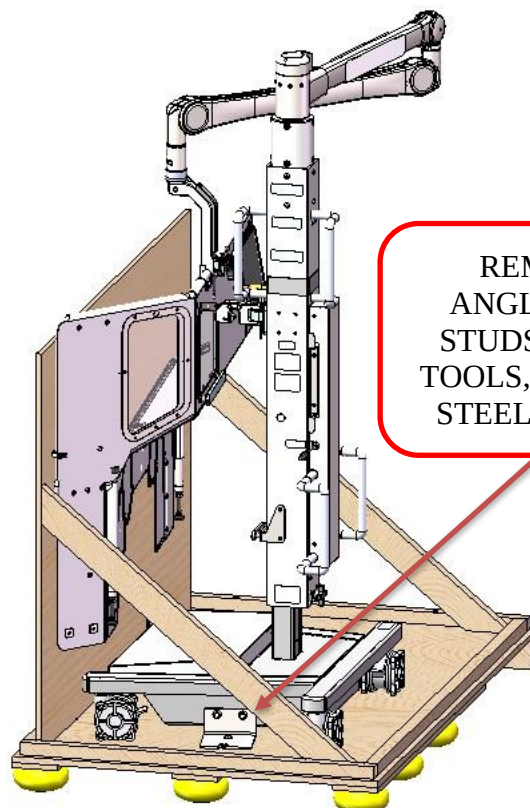


4.


5. Make sure Barrier Pipe is in White Park Block with black Velcro strap around extension before moving off pallet skid. The Barrier bag and green bubble wrap can stay on barrier – IDI Rep. will remove.

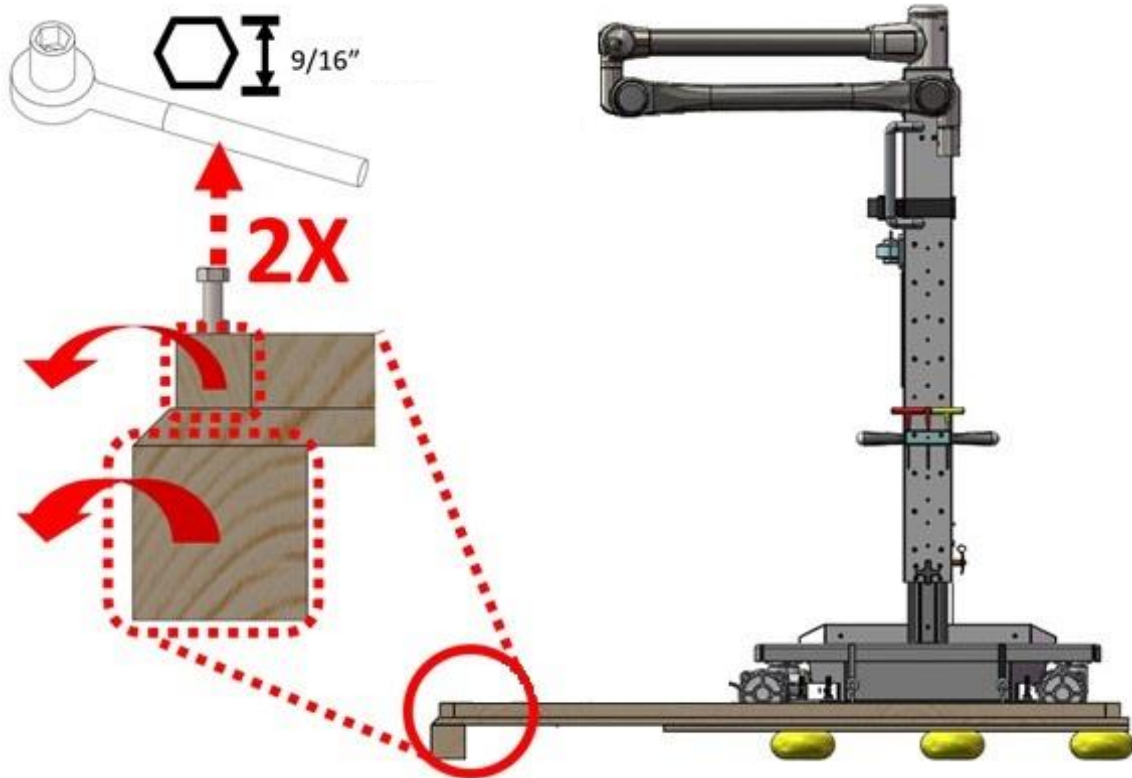

6.

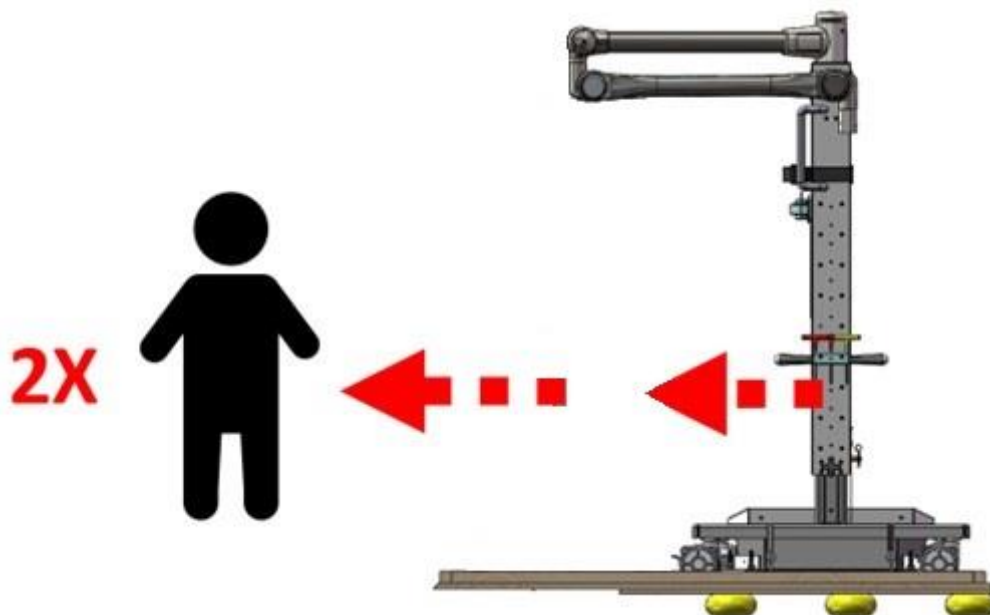

7.

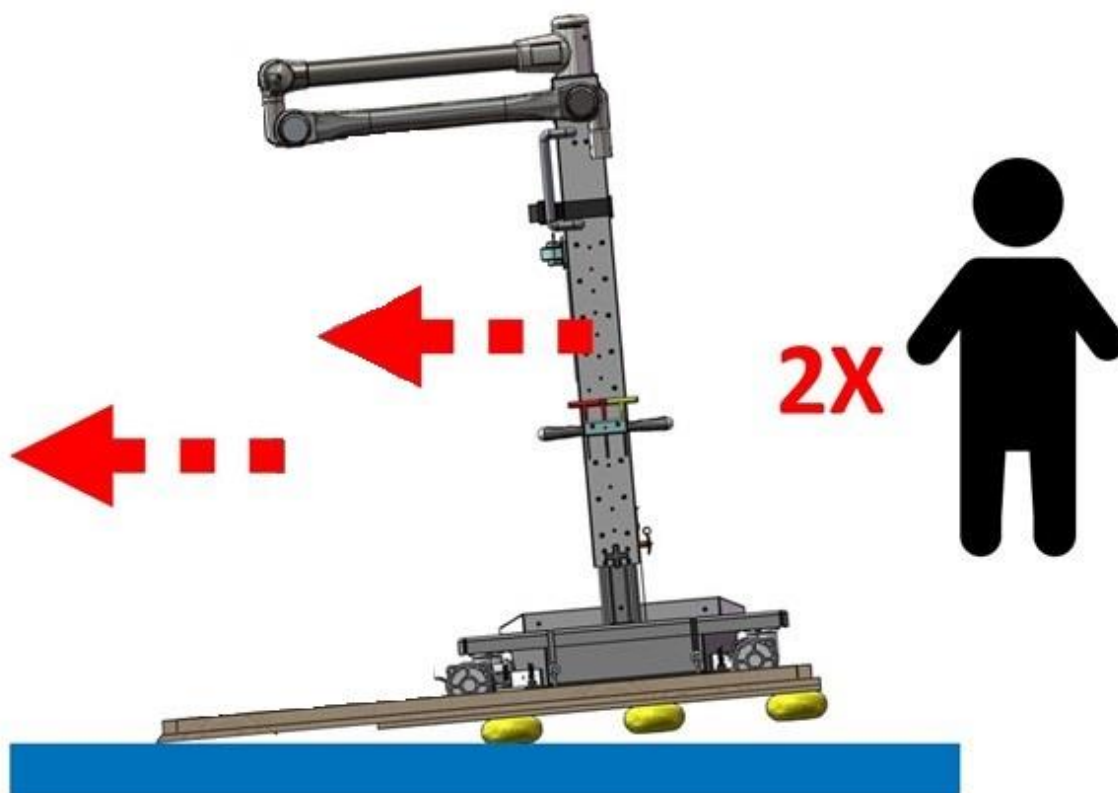

8.

9. After Mobile Stand is off pallet, set the castor lock on for the green labeled castor to lock position to steer to designated location where boxes will be.



10. Remove crating material – Pallet, wood, packing material.

11. Uncrating completed.