


SIMEON HighLine Excellent Light Performance



my.GRIP



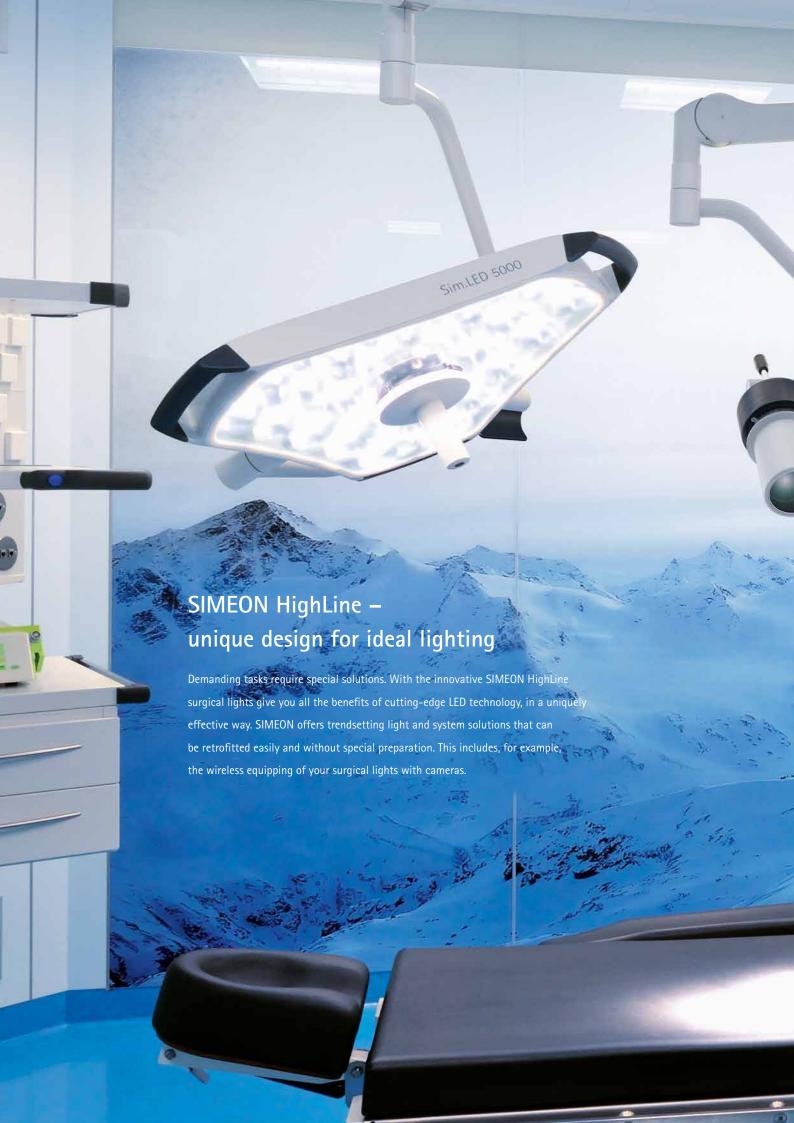



SIMEON HighLine Sim.LED 7000, 5000 – clear advantages

The newest multi-color or single-color LEDs, the patented reflector technology, and the arrangement of the LEDs in a unique design provide more light without shadows. HighLine lights provide optimal illumination for conducting operations.

High-end materials:

- Use of aluminum and safety glass
- Aluminum light head as ideal heat sinks
- Optimal heat management
- Ideal laminar air flow conditions
- Intelligent LED technology: patented Sim.POD reflector

Wireless solution

- Wireless Sim.CAM HD and 4K*, easy replacement and retrofitting of existing surgical lights
- Sim.LED mobile lights with removable battery

Cost effectiveness

- Very long life time > 60,000 h
- Exceptional energy efficiency, low power consumption
- Optimal heat management

Sim.LED 5000

Highest level of hygiene

Innovative antibacterial coating of the light head

Anti.BAC°

- Better hygienic conditions
- Safety glass for ideal hygiene

NEW:

 my.GRIP: sterile control of two operating functions in parallel and directly on the sterile handle

Intelligent solutions

- Innovative Sim.FLEX suspension system
- All surgical lights ready for Sim.CAM HD and Sim.CAM 4K
- my.LED: personalized settings can be simply saved and called up at the press of a button

Cutting-edge LED technology. Revolutionary design. Absolutely high-end.

Patented reflector technology

- Sim.POD: three high-performance LEDs of the latest generation, positioned in the center of a high-gloss aluminum reflector
- sun.LED: The combination of LED and reflector achieves a constantly high color rendering index Ra > 96 in every color temperature setting, with 3,500 K, 4,000 K, 4,500 K, 5,000 K as well as 5,500 K.

Shadow-free light

 Reflector Shadow Resolution System (RSRS) for a unique light quality: large light-emitting surface, shadow-free light, for ideal illumination entirely without sensor systems or automatic adaptation of the light values

Long life time

- > 60,000 h
- 5-year warranty for Sim.PODs
- Guaranteed life time thanks to ideal heat management



- The unique SIMEON triangle design offers a maximum light-emitting surface and an optimal depth of illumination.
- The light body and suspension system, made of light aluminum, allow for a fast and smooth positioning.
- Thanks to the unique light head design,
 allows you to position the lights directly over
 the surgical field virtually without a gap.

User-friendly

 my.LED: personalized settings can be simply saved and called up at the press of a button

NEW:

 my.GRIP: sterile control of two operating functions in parallel and directly on the sterile handle

5-year warranty for Sim.PODs

Environmentally friendly

- Highest energy efficiency
- Minimal heat build-up

SIMEON HighLine – better hygienic conditions and user-friendly operation

Better hygiene thanks to

Anti.BAC®

innovative antibacterial coating

- Highest level of antimicrobial efficacy
 (log 5.9 for Escherichia coli and Staphylococcus aureus)
- Reduces germs by > 99.999 %
- No potentially physiologically harmful nano-silver or free silver ions
- Antimicrobial effect in accordance with ISO 22196 (JIS Z 2801)

Better hygiene thanks to

high-end materials - aluminum and safety glass

- Aluminum light head as the ideal heat sink
- Scratch-resistant safety glass guarantees an optimal light transmittance throughout the entire lifecycle of the lights – no plastic effect!

Sim.LED 7000 MC



NEW:


User-friendly thanks to my.GRIP - One Handle. Personalized Options. SIMEON Multifunction Handle

- Sterile control of two operating functions in parallel and directly on the sterile handle
- Up to 12 freely selectable combinations
- Additional step for optimal hygienic conditions
- Easily programmable using cardanic control panel

- The light intensity, light field and color temperature can all be easily set using the cardanic or wall control
- For endoscopic interventions, the Sim.BIANCE mode is available
- my.LED: personalized settings can be simply saved and called up at the press
 of a button




SIMEON HighLine surgical lights – ideal for the OR and special high-end applications

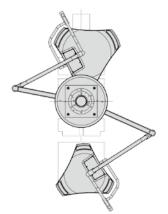


General surgery

A combination of the HighLine Sim.LED 7000/5000 is ideally suited for this application.

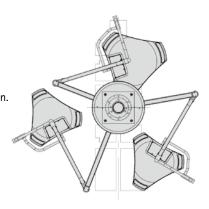

Cardiac OR

A triple combination of HighLine Sim.LED 7000/7000/5000 brings the perfect light for any cardiological intervention.

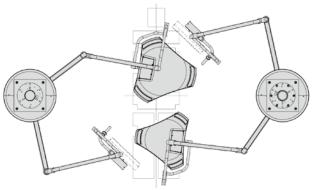

Hybrid OR

A decentralized OR situation will shed light on any situation. This allows other equipment to be used in the middle of the OR table.

The unique design allows for the ideal combination of a wide variety of Sim.LEDs.




Sim.LED 7000/5000

Sim.LED 7000/7000





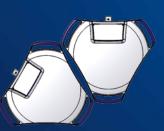


SIMEON HighLine	Sim.LED 7000 SC	Sim.LED 5000 SC	Sim.LED 7000 MC	Sim.LED 5000 MC
Central illumination intensity at a distance of 1 meter [lx]	160,000	140,000	160,000	140,000
Electronical dimmability from/to [lx]	48,000 – 160,000	42,000 – 140,000	48,000 – 160,000	42,000 – 140,000
Color temperature [K]	4,500	4,500	3,500/4,000/ 4,500/5,000/ 5,500	3,500/4,000/ 4,500/5,000/ 5,500
Focusable field size d10 at a distance of 1 meter [mm]	180 – 300	180 – 280	180 – 300	180 – 280
Electronic field adjustment	yes	yes	yes	yes
Color rendering index Ra in every color temperature setting	96	96	96	96
Red rendering index R9	96	96	94	94
"Skin color" rendering index R13	98	98	98	98
Illumination depth without refocusing (L1+L2) at 20% [mm]	1,100	1,200	1,100	1,200
Illumination depth without refocusing (L1+L2) at 60% [mm]	500	600	500	600
Light head weight [kg]	18	15	18	15
Light head power consumption [W]	66	66	60	50
LED life time [h]	> 60,000	> 60,000	> 60,000	> 60,000
Temperature increase at head height	< 1 °C	< 1 °C	< 1 °C	< 1 °C
Light-emitting surface	2,940 cm ²	2,040 cm ²	2,940 cm ²	2,040 cm ²
Sim.CAM	optional	optional	optional	optional
Approval	CE	CE	CE	CE







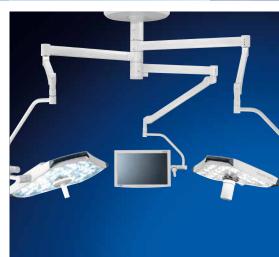

Sim.FLEX

Maximum flexibility for OR planners: With an extension arm length of up to 2,000 mm, our light systems are also easy to use in hybrid ORs. And with Sim.FLEX, it is possible to expand the system by adding mounting arms at any time – even after installation.

Standard suspension system

Easy positioning of the products is an important advantage of the suspension system. Using the standard suspension system, this is possible in single up to quadruple combinations.









Service

Easily accessible and easy to reach electronics.

Controls

SIMEON HighLine OR lights can be simultaneously controlled via wall control or directly at the cardanic control.

Sim.LED mobile lights

- Mobile and flexible:OR lights with an excellent performance
- Mobile stand with one or two batteries
- Ready for use 24/7
- The only mobile wireless solution on the market




NEW: Sim.CAM 4K Wireless* – Wireless* – More Details. Better Treatment.

Fascinating images—with clear, sharp contours and strikingly improved color. Details that Full HD was unable to depict are now clearly and distinctly visualized in 4K ULTRA HD. If needed, a cable version of Sim.CAM 4K as a separate arm solution is available.

Sim.CAM 4K and Sim.CAM 4K Wireless*		
Resolution	8,847,360 pixels	
4K Ultra HD	4096 x 2160	
Zoom	12x optical 12x digital	
Image positioning	> 360°, mechanically and electronically	
Video outputs	HDMI 1.4, DVI	

Technical specifications are subject to change

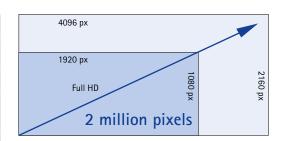

Use 4K ULTRA HD in surgical settings for detail-optimized imagery. The high image quality offers you finely structured images with amazing details and thus a brilliant visual representation.

^{*} Available starting 3rd quarter/2018



More Details. Great Flexibility.

With the new Sim.CAM HD Wireless you obtain a premium HD camera that sends images to your monitor without additional wiring, thanks to the new wireless technology. The integrated HD camera transfers all data sharply and precisely from every selected position. If needed, it is also available as a cable-based solution.

Sim.CAM HD and Sim.CAM HD Wireless		
Resolution	2,000,000 pixels	
Full HD	1,920 x 1,080	
Aspect ratio	16:9	
Zoom	10x optical 12x digital	
Image positioning	> 360°, mechanically and electronically	
Video outputs	HDMI, DVI	

Technical specifications are subject to change





Sim.SCREEN

Sim.SCREEN is also available as a double monitor holder.



NEW:

Sim.Screen 32" monitor holder

Despite the fact that commercial monitors are larger and heavier, the mount's new design guarantees a simple and easy use.







NDS Expand OR

- Medically-approved full HD video recorder and streamer
- Bidirectional audio and video streaming
- Simultaneous control of the NDS Expand OR from the device or through the PC software

Sim.INTERFACE

With the help of the Sim.INTERFACE, all integration systems that are available on the market can control our products.







Excellence in medical care

For over 15 years now, LED technology pioneer SIMEON has provided complete, future-oriented top products at a cost-to-benefit ratio which is unparalleled across the world. Thanks to their special design and unique reflector technology, SIMEON's innovative solutions ensure ideal lighting conditions and cover all demands for hygiene and function in the medical field. SIMEON products are made in Germany and undergo certified production with continuous quality control.

Not all accessories/options shown are available in the U.S.A.



S.I.M.E.O.N. Medical GmbH & Co. KG In Grubenäcker 18 78532 Tuttlingen, Germany Telephone: +49 7461 90068-0

Fax: +49 7461 90068-900

www.simeonmedical.com info@simeonmedical.com

