



# IDI Product Catalog

*Intuitive by design*



idi

Image Diagnostics Inc.



Veteran  
Owned

# Table of Contents

## Imaging Tables and options

Aspect ISR G3 Table	3-7
Aspect 100-4T G3 Table	8-12
Aspect 100-4 Table	13-15
Aspect 100RTL Table	16-19
Aspect 100RT Table	20-23
Aspect 100T Table	24-27
Aspect 100E Table	28-31
Aspect 100UC <i>Plus</i> Urology Table	32-36
Aspect 100US Gyn Table	37-38

## Video Integration Systems and options

MDS	39-42
ilex55	43-44
ilex 32 4K	45
ilex 32	46
ilex C	47
Connexus Direct	48

## SIMEON O.R. Lighting Systems

SIMEON Sim. LED 4500	49
SIMEON Highline & Businessline O.R. Lighting Systems	50

## Comparison Matrices

Aspect Four-Way Float Top Tables Specifications & Accessory Matrix	51
Aspect Urology & GYN Table Specifications & Accessory Matrix	52
Aspect Pain Management Table Specifications & Accessory Matrix	53

## Customer Service

IDI Customer Service Contacts	Back page
-------------------------------	-----------

For product warranties, please request document D000-2069 from an IDI representative.  
Extended warranties are available.

## Aspect ISR G3 Table



### Aspect ISR G3 Table Description:

Mobile Imaging Table with Motorized Elevation, Trendelenburg Tilt, Iso Roll and Manual 4-way Table-top Float. For Advanced Vascular Surgery and Interventional Procedures.

### Aspect ISR G3 #A100-2603 (U.S. accessory rails ) #A100-2834 (Metric accessory rails)

For Advanced Vascular Procedures; Includes short, recessed side rails 5" (127mm) in length at imaging end of table that facilitate use of removable headrest extension (option) and also facilitates attachment of catheter tray extension at narrow end of table. Two separate catheter tray extensions are available: One for the pedestal end of table and one for the imaging end of table. Up to 5 day battery backup pack, 600 lb. patient load, steer assist caster, 3.75" tongue height, 32" longitudinal tabletop float, and 33" to 45" tabletop elevation.

Includes the following:

- A100-2170 Pendant hand control
- X100-2492: Tabletop 22"/559mm x 24"/610mm X 84"/2134mm
- X100-2310 Tabletop Pad, 2" (51mm) thick
- A270-0163: Panning hand control
- C000-0328: Patient restraint straps



## Aspect ISR G3 Table preferred accessory packages

### Aspect ISR G3 Practice Starter Kit

(Peripherals - head/imaging end tray) # A100-2700    Optional package includes:

1. A100-2322 Catheter Tray Tabletop Extension – **Mounts to head/imaging end of table**
2. A100-1007 Clamp-on 8-1/2" US Accessory Rail
3. A310-056 Armboard, quick release, rail mounted, includes pad and strap
4. A100-2244 VAB Carbon Fiber Arm Board w/ pad (Fistula Applications)
5. C000-1181 Medical Grade Anti-Fatigue Mat 24" x 66" x 3/4" thick for maximum comfort & support
6. A610-051 Radiation Shield, Articulating, rail mounted
7. C000-1324 IV Pole, Accessory Rail Mounted, Stainless Steel, 4 bottle hooks
8. C000-1325 Rail clamp, Accessory, Aluminum, for IV Pole

### Aspect ISR G3 Practice Starter Kit - Cardiology # A100-2975    Optional package includes:

1. A100-1769 Catheter Tray Tabletop Extension – **Mounts to Foot/pedestal end of table**
2. A100-1007 Clamp-on 8-1/2 US Accessory Rail
3. A310-056 Armboard: quick-release, rail mounted, includes pad and strap
4. A100-2244 VAB Carbon Fiber Arm Board w/ pad (Fistula Applications)
5. C000-1181 Medical Anti-Fatigue 24" X 66" X 3/4" thick for maximum comfort & support
6. A610-051 Radiation Shield, articulating, rail-mounted
7. C000-1324 IV Pole, Accessory Rail Mounted, Stainless Steel, 4 bottle hooks
8. C000-1325 Rail Clamp, Accessory Aluminum for IV Pole

### Aspect ISR G3 Accessory Kit - Vascular # A100-2977    Optional package includes:

1. A100-2322 Catheter Tray Tabletop Extension – **Mounts to head/imaging end of table**
2. A100-1007 Clamp-on 8-1/2" US accessory rail
3. A310-056 Arm board: quick-release, rail mounted, includes pad and strap
4. A100-2244 VAB Carbon Fiber Arm Board w/ pad (Fistula Applications)









### Aspect ISR G3 Accessory Kit - Cardiology # A100-2978    Optional package includes:

1. A100-1769 Catheter Tray Tabletop Extension – **Mounts to foot/pedestal end of table**
2. A310-056 Armboard – Rail Mount
3. A100-1007 Clamp-on 8-1/2" US Accessory Rail
4. A100-2244 VAB Carbon Fiber Arm Board w/ pad (Fistula Applications)

## Aspect ISR G3 Table Optional and Replacement Accessories:

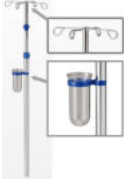

Photo	Description <b>Red text = Metric accessory rails</b>	Catalog Number
	<b>Replacement Tabletop Pad</b> 2" thick foam pad. 19.5"/22"/24" wide x 84" long	X100-2310
	<b>Clamp-on Accessory Rail</b> Rail is 8 ½" long x ¾" thick x 1-1/8" high. <b>Clamp-on Accessory Rail</b> 10mm x 25mm x 216mm clamp-on accessory rail.	A100-1007 A100-2379
	<b>Armboard</b> Quick-release, rail mounted (requires Clamp-on Accessory Rail A100-1007 or A100-2379 above) unless used on a table with built-in accessory rails.	A310-056 A310-062
	<b>Replacement pad for A310-056 or A310-062 Armboard</b>	X100-1755
	<b>VAB Carbon Arm Board w/ Pad</b> Carbon fiber, one-piece, extra wide armboard for vascular access and Fistula application procedures (body section slides under patient/table pad). Includes strap and pad.	A100-2244
	<b>Replacement Pad for A100-2244 VAB Armboard</b>	X100-2243
	<b>Armboard</b> Carbon fiber, shoulder mount (attaches by sliding thin section under patient/table pad), flat. Includes pad and strap.	A310-059
	<b>Replacement Pad for A310-059 Armboard</b>	C000-0689
	<b>Headrest Extension</b> Articulating carbon fiber headrest. Attaches to head end of ISR G3 table. Includes built in Anesthesia Screen Holder sockets. 23-½"/ 587mm wide x 12"/305mm long.	A100-2317 A100-2505
	<b>Replacement Pad for A100-2317 or A100-2505 Headrest Extension</b>	X100-2325
	<b>Peripheral Extension</b> Carbon Fiber 12"/30.5cm table extension for imaging end of table. Table top must have mini rails as shown.	A100-2486
	<b>Replacement Pad for A100-2486 Peripheral Extension</b>	X100-2503
	<b>Tabletop Catheter Tray Extension</b> 24.25"/616mm wide, 30"/762mm long. <b>Mounts to pedestal end of table only.</b> The maximum weight this extension can support is 20 LBS or 9 Kgs.	A100-1769

## Aspect ISR G3 Table Optional and Replacement Accessories:

Photo	Description <span style="color: red;">Red text = Metric accessory rails</span>	Catalog Number
	<b>Tabletop Catheter Tray Extension</b> 20 3/4"/527mm x 27"/686mm long. <b>Mounts to imaging end of table only.</b> Requires tabletop X100-2492. The maximum weight this extension can support is 20 LBS or 9 Kgs.	A100-2322
	<b>Radiation Shield - .5mm lead equivalence</b> Articulating, rail-mounted. Attaches to table side-rail at base end of table, either side, adjustable. Articulating joint allows shield to conform to table contour. Snap-on patient side upper section is removable. 41-3/8" wide x 32-1/8" high x 2-5/8" deep.	A610-051 <span style="color: red;">A610-106</span>
	<b>Anesthesia Screen Holder</b> Flexible "hoop" style. Mounts directly to A100-2317 Headrest Extension. For use on accessory rails - a pair of rotating side rail clamps, part number C000-0746, are required. Rotating side rail clamps and clamp-on accessory rails are not required if table is ordered with the A100-2317 Headrest Extension.	C000-0965
	<b>Pendant Hand Control</b> For Aspect ISR. Controls elevate and tilt functions only (Included with table).	A100-2170
	<b>Panning Hand Control</b> For Aspect ISR G3. (Included with table).	A270-0163 <span style="color: red;">A270-0178</span>
	<b>ISR Table Foot Switch with Safety Hood.</b> Foot control for Aspect ISR G3 table. Foot control motions are: iso-roll, transverse travel, elevate/lower and Trendelenburg down/up. (Safety hood helps avoid accidental motorized table movement from foot switch).	A100-2071
	<b>Heavy Duty Medical Grade Anti-Fatigue Mat</b> (24" (610mm) x 66" (1676mm) x 3/4" (19mm) IDI heavy-duty anti-fatigue mat. One-piece construction for high-traffic areas. Safe & non-toxic, contains no plasticizers and will not leach over time. Durable brushed texture.	C000-1181
	<b>Side-rail clamp</b> Rotating, with socket for attachment of accessories with 5/8 inch (15.88mm) circular post. Requires A100-1007 (or A100-2379 for metric rails) clamp-on accessory rail.	C000-0746 <span style="color: red;">C000-1065</span>



## Aspect ISR G3 Table Optional and Replacement Accessories:

Photo	Description <span style="color: red;">Red text = Metric accessory rails</span>	Catalog Number
	<b>IV Pole</b> Stainless steel with four bottle hooks Range: 1070 -1660 mm Tube diameter: 25 / 18 mm  <b>IV Pole Clamp C000-1325 is required</b>	C000-1324
	<b>IV Pole Clamp</b> Aluminum. 25 mm pole opening	C000-1325

## Shipping Crate Specifications

- Packaged weight (approx.): 890 lbs. (404kg)
- Length: 96" (2438mm).
- Width: 42" (1067mm).
- Height: 47" (1194mm).

## Aspect 100-4T G3 Table



### Aspect 100-4T G3 Table Description:

**Mobile Imaging Table with Motorized Elevation, Trendelenburg Tilt and manual 4-way tabletop float.**

Aspect 100-4T G3 # A100-2602 (U.S accessory rails)

# A100-2831 (Metric accessory rails)

For Advanced Vascular Procedures; Includes short, recessed side rails 5" (127mm) in length at imaging end of table that facilitate use of removable headrest extension (option) and also facilitates attachment of catheter tray extension at narrow end of table. Two separate catheter tray extensions are available: One for the pedestal end of table and one for the imaging end of table. Up to 5 day battery backup pack, 600 lb. patient load, steer assist caster, 3.75" tongue height, 32" longitudinal tabletop float and 32" to 45" tabletop elevation.

Includes the following:

- A100-2190 pendant hand control
- X100-2492: Tabletop 22"/559mm x 24"/610mm X 84"/2134mm
- X100-2310 Tabletop Pad, 2" (51mm) thick.
- A270-0173: Panning hand control
- C000-0328: Patient restraint straps





## Aspect 100-4T G3 Table preferred accessory packages

### Aspect 100-4T G3 Practice Starter Kit

(Peripherals - head/imaging end tray) # A100-2700 Optional package includes:

1. A100-2322 Catheter Tray Tabletop Extension – **Mounts to head/imaging end of table**
2. A100-1007 Clamp-on 8-1/2" US Accessory Rail
3. A310-056 Armboard, quick release, rail mounted, includes pad and strap
4. A100-2244 VAB Carbon Fiber Arm Board w/ pad (Fistula Applications)
5. C000-1181 Medical Grade Anti-Fatigue Mat 24" x 66" x 3/4" thick for maximum comfort & support
6. A610-051 Radiation Shield, Articulating, rail mounted
7. C000-1324 IV Pole, Accessory Rail Mounted, Stainless Steel, 4 bottle hooks
8. C000-1325 Rail clamp, Accessory, Aluminum, for IV Pole

### Aspect 100-4T G3 Practice Starter Kit - Cardiology # A100-2975 Optional package includes:

1. A100-1769 Catheter Tray Tabletop Extension - **Mounts to Foot/pedestal end of table**
2. A100-1007 Clamp-on 8-1/2 US Accessory Rail
3. A310-056 Armboard: quick-release, rail mounted, includes pad and strap
4. A100-2244 VAB Carbon Fiber Arm Board w/ pad (Fistula Applications)
5. C000-1181 Medical Anti-Fatigue 24" X 66" X 3/4" thick for maximum comfort & support
6. A610-051 Radiation Shield, articulating, rail-mounted
7. C000-1324 IV Pole, Accessory Rail Mounted, Stainless Steel, 4 bottle hooks
8. C000-1325 Rail Clamp, Accessory Aluminum for IV Pole

### Aspect 100-4T G3 Accessory Kit - Vascular # A100-2977 Optional package includes:

1. A100-2322 Catheter Tray Tabletop Extension – **Mounts to head/imaging end of table**
2. A100-1007 Clamp-on 8-1/2" US accessory rail
3. A310-056 Arm board: quick-release, rail mounted, includes pad and strap
4. A100-2244 VAB Carbon Fiber Arm Board w/ pad (Fistula Applications)









### Aspect 100-4T G3 Accessory Kit - Cardiology # A100-2978 Optional package includes:

1. A100-1769 Catheter Tray Tabletop Extension – **Mounts to foot/pedestal end of table**
2. A310-056 Armboard – Rail Mount
3. A100-1007 Clamp-on 8-1/2" US Accessory Rail
4. A100-2244 VAB Carbon Fiber Arm Board w/ pad (Fistula Applications)



## Aspect 100-4T G3 Table Optional and Replacement Accessories:

Photo	Description <span style="color: red;">Red text = Metric accessory rails</span>	Catalog Number
	<b>Replacement Tabletop Pad</b> 2" thick foam pad. 19.5"/22"/24" wide x 84" long	X100-2310
	<b>Clamp-on Accessory Rail</b> Rail is 8 1/2" long x 3/8" thick x 1-1/8" high. <span style="color: red;"><b>Clamp-on Accessory Rail</b> 10mm x 25mm x 216mm clamp-on accessory rail.</span>	A100-1007 <span style="color: red;">A100-2379</span>
	<b>Armboard</b> Quick-release, rail mounted (requires Clamp-on Accessory Rail A100-1007 or A100-2379 above) unless used on a table with built-in accessory rails.	A310-056 <span style="color: red;">A310-062</span>
	<b>Replacement pad for A310-056 or A310-062 Armboard</b>	X100-1755
	<b>VAB Carbon Armboard w/ Pad</b> Carbon fiber, one-piece, extra wide armboard for vascular access procedures (body section slides under patient/table pad). Includes strap and pad.	A100-2244
	<b>Replacement Pad for A100-2244 VAB Armboard</b>	X100-2243
	<b>Armboard</b> Carbon fiber, shoulder mount (attaches by sliding thin section under patient/table pad), flat. Includes pad and strap.	A310-059
	<b>Replacement Pad for A310-059 Armboard</b>	C000-0689
	<b>Headrest Extension</b> Articulating carbon fiber headrest. Attaches to head end of ISR G3 table. Includes built in Anesthesia Screen Holder sockets. 23-1/8" / 587mm wide x 12"/305mm long.	A100-2317 <span style="color: red;">A100-2505</span>
	<b>Replacement Pad for A100-2317 or A100-2505 Headrest Extension</b>	X100-2325
	<b>Peripheral Extension</b> Carbon Fiber 12"/30.5cm table extension for imaging end of table. Table top must have mini rails as shown.	A100-2486
	<b>Replacement Pad for A100-2486 Peripheral Extension</b>	X100-2503
	<b>Tabletop Catheter Tray Extension</b> 24.25"/616mm wide, 30"/762mm long. <b>Mounts to pedestal end of table only.</b> The maximum weight this extension can support is 20 LBS or 9 Kgs.	A100-1769

## Aspect 100-4T G3 Table Optional and Replacement Accessories:

Photo	Description <span style="color: red;">Red text = Metric accessory rails</span>	Catalog Number
	<b>Tabletop Catheter Tray Extension</b> 20 3/4"/527mm x 27"/686mm long. <b>Mounts to imaging end of table only.</b> Requires tabletop X100-2492. The maximum weight this extension can support is 20 LBS or 9 Kgs.	A100-2322
	<b>Radiation Shield - .5mm lead equivalence</b> Articulating, rail-mounted. Attaches to table side-rail at base end of table, either side, adjustable. Articulating joint allows shield to conform to table contour. Snap-on patient side upper section is removable. 41-3/8" wide x 32-1/8" high x 2-5/8" deep.	A610-051 <span style="color: red;">A610-106</span>
	<b>Anesthesia Screen Holder</b> Flexible "hoop" style. Mounts directly to A100-2317 Headrest Extension. For use on accessory rails - a pair of rotating side rail clamps, part number C000-0746, are required. Rotating side rail clamps and clamp-on accessory rails are not required if table is ordered with the A100-2317 Headrest Extension.	C000-0965
	<b>Pendant Hand Control</b> For Aspect 100-4T G3. Controls elevate and tilt functions only (Included with table).	A100-2190
	<b>Panning Hand Control</b> For Aspect 100-4T G3. (Included with table).	A270-0173 <span style="color: red;">A270-0179</span>
	<b>100-4T G3 Table Foot Switch with Safety Hood.</b> Foot control for Aspect 100-4T G3 table. Foot control motions are: transverse travel, elevate/lower and Trendelenburg down/up. (Safety hood helps avoid accidental motorized table movement from foot switch).	A100-2360
	<b>Heavy Duty Medical Grade Anti-Fatigue Mat</b> (24"(610mm) x 66"(1676mm) x 3/4" (19mm) IDI heavy-duty anti-fatigue mat. One-piece construction for high-traffic areas. Safe & non-toxic, contains no plasticizers and will not leach over time. Durable brushed texture.	C000-1181
	<b>Side-rail clamp</b> Rotating, with socket for attachment of accessories with 5/8 inch (15.88mm) circular post. Requires A100-1007 (or A100-2379 for metric rails) clamp-on accessory rail.	C000-0746 <span style="color: red;">C000-1065</span>


## Aspect 100-4T G3 Table Optional and Replacement Accessories:

Photo	Description <span style="color: red;">Red text = Metric accessory rails</span>	Catalog Number
	<b>IV Pole</b> Stainless steel with four bottle hooks Range: 1070 -1660 mm Tube diameter: 25 / 18 mm  <b>IV Pole Clamp C000-1325 is required</b>	C000-1324
	<b>IV Pole Clamp</b> Aluminum. 25 mm pole opening	C000-1325

## Shipping Crate Specifications

- Packaged weight (approx.): 890 lbs. (404kg)
- Length: 96" (2438mm).
- Width: 42" (1067mm).
- Height: 47" (1194mm).
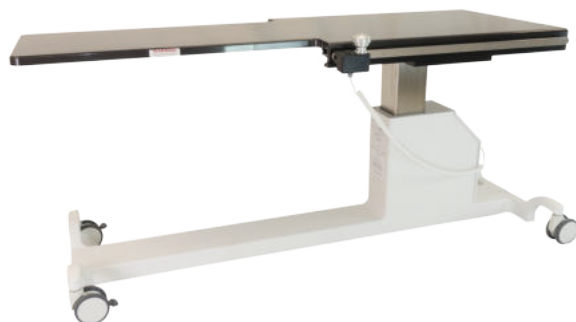
## Aspect 100-4 Table



### Aspect 100-4 Description:

The Aspect 100-4 is designed for C-arm imaging in endovascular or other interventional procedures where a 4-way floating tabletop is required in a non-tilting table.

- 4-Way Manual Floating Tabletop
- Motorized Elevation
- 450 lb. maximum patient weight
- 29.9" to 39.9" vertical travel
- 40" of longitudinal float
- 10" transverse float
- Panning control module
- 120VAC is standard.
- 230VAC is optional and must be specified at time of order.

### Aspect 100-4 #A100-2737 (U.S Rails)    A100-2632 (Metric Rails)

Includes the following:






- X100-2599: Tabletop, 21" to 30" wide (533mm to 762mm) x 80" long (2032mm)
- X100-2650: Tabletop pad, 2" (51mm) thick, for 21" to 30" (533mm to 762mm) x 80" (2032mm) tabletop.
- A270-126 Panning hand control.
- C000-0328: Patient restraint straps.



## Aspect 100-4 Table Optional and Replacement Accessories:

Photo	Description <b>Red text = Metric accessory rails</b>	Catalog Number
	Replacement Tabletop Pad 2" thick foam pad 21"/30" wide x 80" long	X100-2650
	<b>Tabletop Catheter Tray Extension</b> Mounts to pedestal end only.	A100-814
	<b>Armboard</b> Quick-release, rail mounted (requires Clamp-on Accessory Rail A100-1007 or A100-2379 above) unless used on a table with built-in accessory rails.	A310-056 <b>A310-062</b>
	<b>Replacement pad for A310-056 or A310-062 Armboard</b>	X100-1755
	<b>Clamp-on Accessory Rail</b> Rail is 8 1/2" long x 3/8" thick x 1-1/8" high. <b>Clamp-on Accessory Rail</b> 10mm x 25mm x 216mm clamp-on accessory rail.	A100-2218 <b>A100-2838</b>
	<b>Armboard</b> Carbon fiber, shoulder mount (attaches by sliding thin section under patient/table pad), flat. Includes pad and strap.	A310-059
	<b>Replacement Pad for A310-059 Armboard</b>	C000-0689
	<b>VAB Carbon Armboard w/ Pad</b> Carbon fiber, one-piece, extra wide armboard for vascular access procedures (body section slides under patient/table pad). Includes strap and pad.	A100-2244
	<b>Replacement Pad for A100-2244 VAB Armboard</b>	X100-2243
	<b>Tabletop Catheter Tray Extension</b> Mounts to pedestal end only.	A100-814

## Aspect 100-4 Table Optional and Replacement Accessories:

Photo	Description <b>Red text = Metric accessory rails</b>	Catalog Number
	<b>Radiation Shield - .5mm lead equivalence</b> Articulating, rail-mounted. Attaches to table side-rail at base end of table, either side, adjustable. Articulating joint allows shield to conform to table contour. Snap-on patient side upper section is removable. 41- $\frac{3}{8}$ " wide x 32- $\frac{1}{2}$ " high x 2- $\frac{5}{8}$ " deep	A610-051 <b>A610-106</b>
	<b>Panning Hand Control for Aspect 100-4</b> (Included with table)	A270-126 <b>A270-183</b>
	<b>Heavy Duty Medical Grade Anti-Fatigue Mat</b> (24"(610mm) x 66"(1676mm) x $\frac{3}{4}$ " (19mm) IDI heavy-duty anti-fatigue mat. One-piece construction for high-traffic areas. Safe & non-toxic, contains no plasticizers and will not leach over time. Durable brushed texture.	C000-1181
	<b>IV Pole</b> Stainless steel with four bottle hooks Range: 1070 -1660 mm Tube diameter: 25 / 18 mm  <b>IV Pole Clamp C000-1325 is required.</b>	C000-1324
	<b>IV Pole Clamp</b> Aluminum. 25 mm pole opening.	C000-1325

## Aspect 100-4 Shipping Crate Specifications

- Weight (approx.): 635 lbs. (288kg)
- Length: 96" (2438mm)
- Width: 42" (1067mm)
- Height: 47" (1194mm)



## Aspect 100-RTL Table



### Aspect RTL Description:

C-Arm table for pain management and other imaging procedures where motorized elevation, tilt, lateral roll and longitudinal motion are required.

- Motorized elevation
- Trendelenburg tilt
- Lateral roll
- Longitudinal motion
- 120VAC is standard
- 230VAC is optional and must be specified at time of order

### Aspect RTL #A100-1254 (120VAC)    A100-2637 (230VAC)

Includes the following:

- X100-995: General Purpose rectangular tabletop 24" x 80" (610mm) x (2032mm).
- X100-1742: Tabletop pad, 2" (51mm) thick, for 24" (610mm) x 80" (2032mm) top
- K000-0197: Hand control
- C000-0328: Patient restraint straps
- C000-0597: Crescent shaped facial cushion
- C000-0598: Disposable covers for crescent-shaped facial cushion, box of 50

### Aspect RTL #A100-2646 (120VAC)    A100-2647 (230VAC)

Includes the following:

- X100-733: Pain Management tabletop with facial cutout 22" (559mm) x 80" (2032mm).
- X100-1757: Tabletop pad, 2" (51mm) thick, for 22" (559mm) x 80" (2032mm) top with facial cutout
- K000-0197: Hand control
- C000-0328: Patient restraint straps
- C000-0597: Crescent shaped facial cushion
- C000-0598: Disposable covers for crescent-shaped facial cushion, box of 50


## Aspect 100RTL #A100-2672 (120VAC) # A100-2673 (230VAC)

Includes the following:







- **X100-1325: G.I. rectangular tabletop 28" (711mm) x 80" (2032mm).**
- X100-2083: Tabletop pad, 2" (51mm) thick, for 28" (711mm) x 80" (2032mm)
- K000-0197: Hand control
- C000-0328: Patient restraint straps
- C000-0597: Crescent shaped facial cushion
- C000-0598: Disposable covers for crescent-shaped facial cushion, box of 50




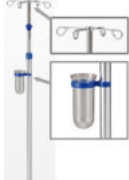

### Aspect 100RTL Table Optional and Replacement Accessories:

Photo	Description <b>Red text = Metric accessory rails</b>	Catalog Number
	<b>Replacement Tabletop Pad</b> For A100-1254 or A100-2637 tables. 2" thick foam pad 24" wide x 80" long.	X100-1742
	<b>Replacement Tabletop Pad</b> For A100-2646 or A100-2647 tables. 2" thick foam pad. 22" wide with a facial cutout x 80" long.	X100-1757
	<b>Replacement Tabletop Pad</b> For A100-2672 or A100-2673 tables. 2" thick foam pad 28" wide x 80" long.	X100-2083
	<b>Side-rail clamp</b> Rotating, with socket for attachment of accessories with 5/8 inch (15.88mm) circular post. Requires A100-1007 (or A100-2379 for metric rails) clamp-on accessory rail.	C000-0746 <b>C000-1065</b>
	<b>Armboard</b> Quick-release, rail mounted (requires Clamp-on Accessory Rail A100-1007 or A100-2379).	A310-056 <b>A310-062</b>
	<b>Replacement pad for A310-056 or A310-062 Armboard</b>	X100-1755
	<b>Clamp-on Accessory Rail</b> Rail is 8 1/2" long x 3/8" thick x 1-1/8" high. <b>Clamp-on Accessory Rail</b> 10mm x 25mm x 216mm clamp-on accessory rail.	A100-1007 <b>A100-2379</b>
	<b>Arm Board</b> Carbon fiber, shoulder mount (attaches by sliding thin section under patient/table pad), flat. Includes pad and strap.	A310-059
	<b>Replacement Pad for A310-059 Armboard</b>	C000-0689
	<b>Hand Control, 100RTL</b> (must verify table serial # with IDI for compatibility). Motions controlled: Elevation, trendelenburg, roll, slide. (Included with table).	K000-0197

## Aspect 100RTL Table Optional and Replacement Accessories:

Photo	Description <span style="color: red;">Red text = Metric accessory rails</span>	Catalog Number
	<b>Foot Control for 100RTL</b> Foot Control unit controls all table motions.	A100-1961
	<b>Scatter Guard radiation shield draping - .5mm lead equivalence</b> For 28" (711mm) wide optional tabletop X100-1325 Dimensions: 59" (1500mm) long, 24" (610mm) below tabletop. Draping extends along both sides and head--end of table. Supported by fabric section over tabletop. <i>Not for use with arm board A310-056.</i>	A610-108
	<b>Radiation Shield - .5mm lead equivalence</b> Fixed shield with integrated clamp system. For use with tables that do not have side rails. Clamps to table sides or at head-end of table. Snap-on upper section above level of tabletop is removable. Use with 100RTL, 100RT, 100T, 100E tables. 24" wide x 32-1/8" high x 4" deep.	A610-080
	<b>Scatter-Guard radiation shield draping - .5mm lead equivalence</b> For for 24" (610mm) rectangular tabletop X100-995. Dimensions: 59" (1500mm) long. 24" (610mm) draping below tabletop. Draping extends along both side and head-end of table. Supported by fabric section over tabletop <i>Note: Not for use with arm board A310-056.</i>	A610-096
	<b>Scatter-Guard radiation shield draping - .5mm lead equivalence</b> for 22" (559mm) tabletop X100-733 models with 22" (559mm) tapered-end tabletop with facial cutout. Dimensions: 59" (1500mm) long. 24" (610mm) draping below tabletop. Draping extends along both side and head-end of table. Supported by fabric section over tabletop. <i>Note: Not for use with arm board A310-056.</i>	A610-092
	<ul style="list-style-type: none"> <li>• Facial Cushion, crescent shape.</li> <li>• Disposable covers for facial cushion, box of 50.</li> </ul>	Cushion: #C000-0597  Covers: #C000-0598

## Aspect 100RTL Table Optional and Replacement Accessories:

	<p><b>Heavy Duty Medical Grade Anti-Fatigue Mat</b> (24"(610mm) x 66"(1676mm) x ¾" (19mm) IDI heavy-duty anti-fatigue mat. One-piece construction for high-traffic areas. Safe &amp; non-toxic, contains no plasticizers and will not leach over time. Durable brushed texture.</p>	<p>C000-1181</p>
	<p><b>IV Pole</b> Stainless steel with four bottle hooks Range: 1070 -1660 mm Tube diameter: 25 / 18 mm <b>IV Pole Clamp C000-1325 is required.</b></p>	<p>C000-1324</p>
	<p><b>IV Pole Clamp</b> Aluminum. 25 mm pole opening  <b>Clamp-on Accessory Rail A100-1007 is required.</b></p>	<p>C000-1325</p>

## 100RTL Shipping Crate Specifications

- Weight (approx.): 605 lbs. (274kg)
- Length: 96" (2438mm).
- Width: 42" (1067mm).
- Height: 47" (1194mm)

## Aspect 100RT Mobile Imaging Table



### Aspect 100RT Description:

C-arm table for pain management and other imaging procedures where motorized elevation, tilt and lateral roll are required.

- Motorized Elevation
- Trendelenburg tilt
- Lateral Roll
- 120VAC is standard
- 230VAC is optional and must be specified at time of order.

### Aspect 100RT #A100-870 (120VAC)    A100-2636 (230VAC)

Includes the following:

- X100-995: General Purpose rectangular tabletop 24" x 80" (610mm) x (2032mm).
- X100-1742: Tabletop pad, 2" (51mm) thick, for 24" (610mm) x 80" (2032mm) top
- K000-0195: Hand control.
- C000-0328: Patient restraint straps.
- C000-0597: Crescent shaped facial cushion.
- C000-0598: Disposable covers for crescent-shaped facial cushion, box of 50.



### Aspect 100RT #A100-2644 (120VAC)    A100-2645 (230VAC)

Includes the following:

- X100-733: Pain Management tabletop with facial cutout 22" (559mm) x 80" (2032mm).
- X100-1757: Tabletop pad, 2" (51mm) thick, for 22" (559mm) x 80" (2032mm) top with facial cutout.
- K000-0195: Hand control.
- C000-0328: Patient restraint straps.
- C000-0597: Crescent shaped facial cushion.
- C000-0598: Disposable covers for crescent-shaped facial cushion, box of 50.










## Aspect 100RT #A100-2660 (120VAC) #A100-2661(230VAC)

Includes the following:


- X100-1325: G.I. rectangular tabletop 28" (711mm) x 80" (2032mm).
- X100-2083: Tabletop pad, 2" (51mm) thick, for 28" (711mm) x 80" (2032mm).
- K000-0195: Hand control.
- C000-0328: Patient restraint straps.
- C000-0597: Crescent shaped facial cushion.
- C000-0598: Disposable covers for crescent-shaped facial cushion, box of 50.



### Aspect 100RT Table Optional and Replacement Accessories:


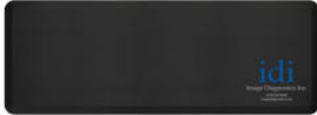

Photo	Description <span style="color: red;">Red text = Metric accessory rails</span>	Catalog Number
	<b>Replacement Tabletop Pad</b> For A100-870 or A100-2636 tables. 2" thick foam pad 24" wide x 80" long.	X100-1742
	<b>Replacement Tabletop Pad</b> For A100-2644 or A100-2645 tables. 2" thick foam pad. 22" wide with a facial cutout x 80" long.	X100-1757
	<b>Replacement Tabletop Pad</b> For A100-2660 or A100-2661 tables. 2" thick foam pad 28" wide x 80" long.	X100-2083
	<b>Side-rail clamp</b> Rotating, with socket for attachment of accessories with 5/8 inch (15.88mm) circular post. Requires A100-1007 (or A100-2379 for metric rails) clamp-on accessory rail.	C000-0746 <span style="color: red;">C000-1065</span>
	<b>Armboard</b> Quick-release, rail mounted (requires Clamp-on Accessory Rail A100-1007 or A100-2379.	A310-056 <span style="color: red;">A310-062</span>
	<b>Replacement pad for A310-056 or A310-062 Armboard</b>	X100-1755
	<b>Clamp-on Accessory Rail</b> Rail is 8 1/2" long x 3/8" thick x 1-1/8" high. <span style="color: red;">Clamp-on Accessory Rail</span> <span style="color: red;">10mm x 25mm x 216mm clamp-on accessory rail.</span>	A100-1007 <span style="color: red;">A100-2379</span>
	<b>Arm Board</b> Carbon fiber, shoulder mount (attaches by sliding thin section under patient/table pad), flat. Includes pad and strap.	A310-059
	<b>Replacement Pad for A310-059 Armboard</b>	C000-0689
	<b>Hand Control, 100RT</b> Must verify table serial # with IDI for compatibility. (Included with table).	K000-0195

## Aspect 100RT Table Optional and Replacement Accessories:

Photo	Description <b>Red text = Metric accessory rails</b>	Catalog Number
	<b>Foot Control</b> for 100RT (must verify table serial # with IDI for compatibility).	A100-1960
	<b>Scatter Guard radiation shield draping - .5mm lead equivalence</b> For 28" (711mm) wide optional tabletop X100-1325 Dimensions: 59" (1500mm) long, 24" (610mm) below tabletop. Draping extends along both sides and head--end of table. Supported by fabric section over tabletop. <i>Not for use with arm board A310-056.</i>	A610-108
	<b>Radiation Shield - .5mm lead equivalence</b> Fixed shield with integrated clamp system. For use with tables that do not have side rails. Clamps to table sides or at head-end of table. Snap-on upper section above level of tabletop is removable. Use with 100RTL, 100RT, 100T, 100E tables. 24" wide x 32-1/8" high x 4" deep.	A610-080
	<b>Scatter-Guard radiation shield draping - .5mm lead equivalence</b> For for 24" (610mm) rectangular tabletop X100-995. Dimensions: 59" (1500mm) long. 24" (610mm) draping below tabletop. Draping extends along both side and head-end of table. Supported by fabric section over tabletop <i>Note: Not for use with arm board A310-056.</i>	A610-096
	<b>Scatter-Guard radiation shield draping - .5mm lead equivalence</b> for 22" (559mm) tabletop X100-733 models with 22" (559mm) tapered-end tabletop with facial cutout. Dimensions: 59" (1500mm) long. 24" (610mm) draping below tabletop. Draping extends along both side and head-end of table. Supported by fabric section over tabletop. <i>Note: Not for use with arm board A310-056.</i>	A610-092
	<ul style="list-style-type: none"> <li>• Facial Cushion, crescent shape.</li> <li>• Disposable covers for facial cushion, box of 50.</li> </ul>	Cushion: #C000-0597 Covers: #C000-0598

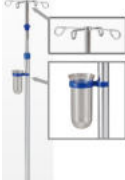

## Aspect 100RT Table Optional and Replacement Accessories:

	<p><b>Heavy Duty Medical Grade Anti-Fatigue Mat</b> (24"(610mm) x 66"(1676mm) x ¾" (19mm) IDI heavy-duty anti-fatigue mat. One-piece construction for high-traffic areas. Safe &amp; non-toxic, contains no plasticizers and will not leach over time. Durable brushed texture.</p>	<p>C000-1181</p>
	<p><b>IV Pole</b> Stainless steel with four bottle hooks Range: 1070 -1660 mm Tube diameter: 25 / 18 mm <b>IV Pole Clamp C000-1325 is required.</b></p>	<p>C000-1324</p>
	<p><b>IV Pole Clamp</b> Aluminum. 25 mm pole opening  <b>Clamp-on Accessory Rail A100-1007 is required.</b></p>	<p>C000-1325</p>

## 100RT Shipping Crate Specifications

- Weight (approx.): 571 lbs. (259 kg.)
- Length: 96" (2438mm)
- Width: 42" (1067mm)
- Height: 47" (1194mm)

## Aspect 100T Mobile Imaging Table



### Aspect 100T Description:

C-arm table for pain management and other procedures where motorized elevation and Trendelenburg tilt are required.

- Motorized Elevation
- Trendelenburg tilt
- 120VAC is standard
- 230VAC is optional and must be specified at time of order.

### Aspect 100T #A100-658 (120VAC)    A100-2638 (230VAC)

Includes the following:

- X100-995: General Purpose rectangular tabletop 24" x 80" (610mm) x (2032mm).
- X100-1742: Tabletop pad, 2" (51mm) thick, for 24" (610mm) x 80" (2032mm) top
- K000-0054: Hand control.
- C000-0328: Patient restraint straps.
- C000-0597: Crescent shaped facial cushion.
- C000-0598: Disposable covers for crescent-shaped facial cushion, box of 50.

### Aspect 100T #A100-2648 (120VAC)    A100-2649 (230VAC)

Includes the following:







- X100-733: Pain Management tabletop with facial cutout 22" (559mm) x 80" (2032mm).
- X100-1757: Tabletop pad, 2" (51mm) thick, for 22" (559mm) x 80" (2032mm) top with facial cutout.
- K000-0054: Hand control.
- C000-0328: Patient restraint straps.
- C000-0597: Crescent shaped facial cushion.
- C000-0598: Disposable covers for crescent-shaped facial cushion, box of 50.

## Aspect 100T #A100-2658 (120VAC) A100-2659 (230VAC)

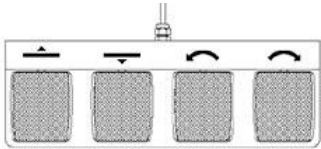





Includes the following:

- X100-1325: G.I. rectangular tabletop 28" (711mm) x 80" (2032mm).
- X100-2083: Tabletop pad, 2" (51mm) thick, for 28" (711mm) x 80" (2032mm).
- K000-0054: Hand control.
- C000-0328: Patient restraint straps.
- C000-0597: Crescent shaped facial cushion.
- C000-0598: Disposable covers for crescent-shaped facial cushion, box of 50.




### Aspect 100T Table Optional and Replacement Accessories:

Photo	Description <span style="color: red;">Red text = Metric accessory rails</span>	Catalog Number
	<b>Replacement Tabletop Pad</b> For A100-658 or A100-2638 tables. 2" thick foam pad 24" wide x 80" long.	X100-1742
	<b>Replacement Tabletop Pad</b> For A100-2648 or A100-2649 tables. 2" thick foam pad. 22" wide with a facial cutout x 80" long.	X100-1757
	<b>Replacement Tabletop Pad</b> For A100-2658 or A100-2659 tables. 2" thick foam pad 28" wide x 80" long.	X100-2083
	<b>Side-rail clamp</b> Rotating, with socket for attachment of accessories with 5/8 inch (15.88mm) circular post. Requires A100-1007 (or A100-2379 for metric rails) clamp-on accessory rail.	C000-0746 <span style="color: red;">C000-1065</span>
	<b>Armboard</b> Quick-release, rail mounted (requires Clamp-on Accessory Rail A100-1007 or A100-2379).	A310-056 <span style="color: red;">A310-062</span>
	<b>Replacement pad for A310-056 or A310-062 Armboard</b>	X100-1755
	<b>Clamp-on Accessory Rail</b> Rail is 8 1/2" long x 3/8" thick x 1-1/8" high. <span style="color: red;">Clamp-on Accessory Rail</span> <span style="color: red;">10mm x 25mm x 216mm clamp-on accessory rail.</span>	A100-1007 <span style="color: red;">A100-2379</span>
	<b>Arm Board</b> Carbon fiber, shoulder mount (attaches by sliding thin section under patient/table pad), flat. Includes pad and strap.	A310-059
	<b>Replacement Pad for A310-059 Armboard</b>	C000-0689
	<b>Hand Control</b> for 100T table. Elevate and tilt functions. (Included with table).	K000-0054

## Aspect 100T Table Optional and Replacement Accessories:

Photo	Description <span style="color: red;">Red text = Metric accessory rails</span>	Catalog Number
	<b>Foot Control for 100T</b> Substitute for hand control. Cannot use hand and foot control at the same time.	A100-500
	<b>Scatter Guard radiation shield draping - .5mm lead equivalence</b> For 28" (711mm) wide optional tabletop X100-1325 Dimensions: 59" (1500mm) long, 24" (610mm) below tabletop. Draping extends along both sides and head--end of table. Supported by fabric section over tabletop. <i>Not for use with arm board A310-056.</i>	A610-108
	<b>Radiation Shield - .5mm lead equivalence</b> Fixed shield with integrated clamp system. For use with tables that do not have side rails. Clamps to table sides or at head-end of table. Snap-on upper section above level of tabletop is removable. Use with 100RTL, 100RT, 100T, 100E tables. 24" wide x 32-1/8" high x 4" deep.	A610-080
	<b>Scatter-Guard radiation shield draping - .5mm lead equivalence</b> For for 24" (610mm) rectangular tabletop X100-995. Dimensions: 59" (1500mm) long. 24" (610mm) draping below tabletop. Draping extends along both side and head-end of table. Supported by fabric section over tabletop <i>Note: Not for use with arm board A310-056.</i>	A610-096
	<b>Scatter-Guard radiation shield draping - .5mm lead equivalence</b> for 22" (559mm) tabletop X100-733 models with 22" (559mm) tapered-end tabletop with facial cutout. Dimensions: 59" (1500mm) long. 24" (610mm) draping below tabletop. Draping extends along both side and head-end of table. Supported by fabric section over tabletop. <i>Note: Not for use with arm board A310-056.</i>	A610-092
	<ul style="list-style-type: none"> <li>• Facial Cushion, crescent shape.</li> <li>• Disposable covers for facial cushion, box of 50.</li> </ul>	Cushion: #C000-0597 Covers: #C000-0598

## Aspect 100T Table Optional and Replacement Accessories:

	<p><b>Heavy Duty Medical Grade Anti-Fatigue Mat</b> (24"(610mm) x 66"(1676mm) x ¾" (19mm) IDI heavy-duty anti-fatigue mat. One-piece construction for high-traffic areas. Safe &amp; non-toxic, contains no plasticizers and will not leach over time. Durable brushed texture.</p>	<p>C000-1181</p>
	<p><b>IV Pole</b> Stainless steel with four bottle hooks Range: 1070 -1660 mm Tube diameter: 25 / 18 mm <b>IV Pole Clamp C000-1325 is required.</b></p>	<p>C000-1324</p>
	<p><b>IV Pole Clamp</b> Aluminum. 25 mm pole opening  <b>Clamp-on Accessory Rail A100-1007 is required.</b></p>	<p>C000-1325</p>

## 100T Shipping Crate Specifications

- Weight (approx): 540 lbs./254kg
- Length: 96" (2438mm)
- Width: 42" (1067mm).
- Height: 47" (1194mm)

## Aspect 100E Mobile Imaging Table



### Aspect 100E Description:

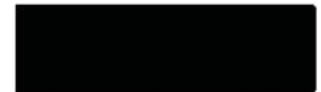
C-arm table for pain management and other procedures where motorized elevation is the only movement required.

- 120VAC is standard.
- 230VAC is optional and must be specified at time of order.

### Aspect 100E #A100-402 (120VAC)    A100-2635 (230VAC)

Includes the following:

- X100-995: General Purpose rectangular tabletop 24" x 80" (610mm) x (2032mm).
- X100-1742: Tabletop pad, 2" (51mm) thick, for 24" (610mm) x 80" (2032mm) top
- K000-0041: Hand control.
- C000-0328: Patient restraint straps.
- C000-0597: Crescent shaped facial cushion.
- C000-0598: Disposable covers for crescent-shaped facial cushion, box of 50.



### Aspect 100E #A100-2642 (120VAC)    A100-2643 (230VAC)

Includes the following:

- X100-733: Pain Management tabletop with facial cutout 22" (559mm) x 80" (2032mm).
- X100-1757: Tabletop pad, 2" (51mm) thick, for 22" (559mm) x 80" (2032mm) top with facial cutout.
- K000-0041 Pendant hand control.
- C000-0328: Patient restraint straps.
- C000-0597: Crescent shaped facial cushion.
- C000-0598: Disposable covers for crescent-shaped facial cushion, box of 50.



## Aspect 100E #A100-2656 (120VAC) A100-2657 (230VAC)

Includes the following:

- X100-1325: GI rectangular tabletop 28" (711mm) x 80" (2032mm).
- X100-2083: Tabletop pad, 2" (51mm) thick, for 28" (711mm) x 80" (2032mm).
- K000-0041: Hand control.
- C000-0328: Patient restraint straps.
- C000-0597: Crescent shaped facial cushion.
- C000-0598: Disposable covers for crescent-shaped facial cushion, box of 50.

### Aspect 100E Table Optional and Replacement Accessories:


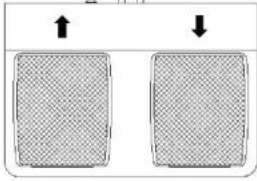






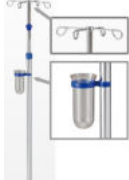

Photo	Description <b>Red text = Metric accessory rails</b>	Catalog Number
	<b>Replacement Tabletop Pad</b> For A100-402 or A100-2635 tables. 2" thick foam pad 24" wide x 80" long.	X100-1742
	<b>Replacement Tabletop Pad</b> For A100-2642 or A100-2643 tables. 2" thick foam pad. 22" wide with a facial cutout x 80" long.	X100-1757
	<b>Replacement Tabletop Pad</b> For A100-2656 or A100-2657 tables. 2" thick foam pad 28" wide x 80" long.	X100-2083
	<b>Side-rail clamp</b> Rotating, with socket for attachment of accessories with 5/8 inch (15.88mm) circular post. Requires A100-1007 (or A100-2379 for metric rails) clamp-on accessory rail.	C000-0746 <b>C000-1065</b>
	<b>Armboard</b> Quick-release, rail mounted (requires Clamp-on Accessory Rail A100-1007 or A100-2379).	A310-056 <b>A310-062</b>
	<b>Replacement pad for A310-056 or A310-062 Armboard</b>	X100-1755
	<b>Clamp-on Accessory Rail</b> Rail is 8 1/2" long x 3/8" thick x 1-1/8" high. <b>Clamp-on Accessory Rail</b> 10mm x 25mm x 216mm <b>clamp-on accessory rail.</b>	A100-1007 <b>A100-2379</b>
	<b>Arm Board</b> Carbon fiber, shoulder mount (attaches by sliding thin section under patient/table pad), flat. Includes pad and strap.	A310-059
	<b>Replacement Pad for A310-059 Armboard</b>	C000-0689
	<b>Hand control for 100E.</b> (Included with table).	K000-0041



Photo	Description <span style="color: red;">Red text = Metric accessory rails</span>	Catalog Number
	<b>Foot Control 100E</b>	A100-758
	<b>Scatter Guard radiation shield draping - .5mm lead equivalence</b> For 28" (711mm) wide optional tabletop X100-1325 Dimensions: 59" (1500mm) long, 24" (610mm) below tabletop. Draping extends along both sides and head--end of table. Supported by fabric section over tabletop. <i>Not for use with arm board A310-056.</i>	A610-108
	<b>Radiation Shield - .5mm lead equivalence</b> Fixed shield with integrated clamp system. For use with tables that do not have side rails. Clamps to table sides or at head-end of table. Snap-on upper section above level of tabletop is removable. Use with 100RTL, 100RT, 100T, 100E tables. 24" wide x 32-1/8" high x 4" deep.	A610-080
	<b>Scatter-Guard radiation shield draping - .5mm lead equivalence</b> For for 24" (610mm) rectangular tabletop X100-995. Dimensions: 59" (1500mm) long. 24" (610mm) draping below tabletop. Draping extends along both side and head-end of table. Supported by fabric section over tabletop <i>Note: Not for use with arm board A310-056.</i>	A610-096
	<b>Scatter-Guard radiation shield draping - .5mm lead equivalence</b> for 22" (559mm) tabletop X100-733 models with 22" (559mm) tapered-end tabletop with facial cutout. Dimensions: 59" (1500mm) long. 24" (610mm) draping below tabletop. Draping extends along both side and head-end of table. Supported by fabric section over tabletop. <i>Note: Not for use with arm board A310-056.</i>	A610-092
	<ul style="list-style-type: none"> <li>• Facial Cushion, crescent shape.</li> <li>• Disposable covers for facial cushion, box of 50.</li> </ul>	Cushion: #C000-0597 Covers: #C000-0598

## Aspect 100E Table Optional and Replacement Accessories:

	<p><b>Heavy Duty Medical Grade Anti-Fatigue Mat</b> (24"(610mm) x 66"(1676mm) x ¾" (19mm) IDI heavy-duty anti-fatigue mat. One-piece construction for high-traffic areas. Safe &amp; non-toxic, contains no plasticizers and will not leach over time. Durable brushed texture.</p>	<p>C000-1181</p>
	<p><b>IV Pole</b> Stainless steel with four bottle hooks Range: 1070 -1660 mm Tube diameter: 25 / 18 mm <b>IV Pole Clamp C000-1325 is required.</b></p>	<p>C000-1324</p>
	<p><b>IV Pole Clamp</b> Aluminum. 25 mm pole opening  <b>Clamp-on Accessory Rail A100-1007 is required.</b></p>	<p>C000-1325</p>

## 100E Shipping Crate Specifications

- Weight (approx.): 495 lbs. (225kg)
- Length: 96" (2438mm)
- Width: 42" (1067mm).
- Height: 47" (1194mm).

## Aspect 100-UCPLUS Urology/Multi-purpose Table



Shown with optional leg holders (#A100-2247).



Shown with tabletop extension (A100-1400)

### Aspect 100UCPlus description:

Imaging table for urology and multi-purpose use. Urological accessories are easily removable so the table can be quickly configured for multi-purpose C-arm imaging or general patient positioning for non-urological procedures. A removable radiographic cassette tray provides KUB radiography compatibility with mobile x-ray machines, including both conventional film/grid systems and digital detectors.




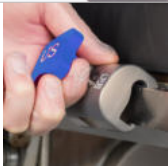

- 120VAC is standard.
- 230VAC is optional and must be specified at time of order.

### Aspect 100UCPlus #A100-1517 (120VAC)    A100-2639 (230VAC)








Includes the following:

- A100-1876 Stirrups Mounting Rail extension Kit
- X100-1818: tabletop pad, 28" (711mm) x 46 ½" (1181mm) x 2" (51mm) thick
- A100-1400 Tabletop extension, 24" (610mm)W x 34" (864mm) (Metric tables #A100-2675)
- X100- 1817: Tabletop extension pad, 24" (610mm) x 34" (864mm) x 2" (51mm) thick
- K000-0158: Pendant hand control
- A100-1540: Foot switch control
- U100-3025: Removable x-ray shield flap, for perineal end of table
- C000-0571: Patient restraint straps
- A310-056: Arm board, rail mounted, with pad and strap (A310-062 for metric rails)
- C000-0595: Urology drain bag support hoop
- C000-0593: Disposable drain bag & hose assembly, sterile, box of 20.
- C000-0492: Disposable covers for foot switch, box of 50
- C000-0645: Disposable tabletop pad protector (qty. of 1)
- C000-0612: 5 gallon collection bag (qty of 1)



## Aspect 100UCPlus Optional and Replacement Accessories:

Photo	Description <b>Red text = Metric accessory rails</b>	Catalog Number
	<b>IDI ComfortFin 350 Lithotomy Leg holders (Stirrups pair only)</b> ComfortFin Lithotomy Leg holders are designed to improve infection control with durable, easy to clean silicone straps. Includes integrated lateral fin to help protect the head of the fibula and the peroneal nerve. 4" more working space with an increased range of motion. Patient weight capacity: 350 lbs (159 kg).	C100-3042
	<b>Replacement Boot with Straps</b> for IDI C100-3042 ComfortFin 350 Lithotomy Leg Holders	C100-3080
	<b>Kit:</b> IDI ComfortFin 350 Lithotomy Leg holders C100-3042 with C100-3046 economy side-rail clamps	A100-2246 <b>A100-3074</b>
	<b>IDI ComfortFin 500 Lithotomy Leg holders (Stirrups pair only)</b> ComfortFin Lithotomy Leg holders are designed to improve infection control with durable, easy to clean silicone straps. Includes integrated lateral fin to help protect the head of the fibula and the peroneal nerve. Boots are sized for today's heavier patients. 4" more working space with an increased range of motion. Patient weight capacity: 500 lbs. (227kg).	C100-3045
	<b>Replacement Boot with Straps</b> for IDI C100-3045 ComfortFin 500 Lithotomy Leg Holders	C100-3079
	<b>Kit:</b> IDI ComfortFin 500 Lithotomy Leg holders C100-3045 with C100-3046 economy side-rail clamps.	A100-2247 <b>A100-3075</b>
	<b>Side-rail clamp, removable. IDI 4100.</b> Attaches and locks anywhere on the side rail for fast and easy mounting. A push button safely releases rail clamp—helps to prevent false rail engagement. Autoclavable.	C100-3043  <b>C100-3076</b>
	<b>Economy side-rail clamp.</b> Attaches anywhere on the table rail for fast and easy mounting. Rounded ergonomic shape is easy on the hands. Autoclavable.	C100-3046  <b>C100-3077</b>
	<b>Lithotomy Leg Holder Cart</b> Easy to push/pull with ergonomic handles. Stirrups are isolated to prevent damage. Clamp storage compartment. Designed to safely store ComfortFin stirrups when not in use.	C100-3044




## Aspect 100UC<sup>Plus</sup> Optional and Replacement Accessories:

Photo	Description <span style="color: red;">Red text = Metric accessory rails</span>	Catalog Number
	<b>Replacement Tabletop Pad</b> 2 inch thick foam pad. 28" wide x 46 1/2" long.	X100-1818
	<b>Replacement Extension Tabletop Pad</b> 2 inch thick foam pad. 28" wide x 34" long.	X100-1817
	<b>Single Use Stirrup Boot Cover</b> Waterproof single use stirrup boot covers help eliminate strike through and bacterial migration. Infectious patients (known or suspected). Patients with open sores, abrasion or lesions. Protect against cross contamination.	C100-3070
	<b>Physician Elbow Supports (pair)</b> Removable	C000-0594
	<b>Armboard</b> Quick-release, rail mounted (requires Clamp-on Accessory Rail A100-1007 or A100-2379 above) unless used on a table with built-in accessory rails.	A310-056 <span style="color: red;">A310-062</span>
	<b>Replacement pad for A310-056 or A310-062 Armboard</b>	X100-1755
	<b>Split Leg Holders</b> For use with Aspect 100UC Plus Tables. Includes pads. Prone percutaneous, lateral or modified lateral nephrostomy. Includes two split leg holder clamps # C000-1115 or <span style="color: red;">C000-1217</span> .	C000-1146 <span style="color: red;">C000-1216</span>
	<b>Right replacement pad for Split Leg Holder</b> C000-1146 or A310-062	C000-1357
	<b>Left replacement pad for Split Leg Holder</b> C000-1146 or A310-062	C000-1358
	<b>Replacement Split Leg Holder Clamp</b> For use with split leg holders. # C000-1146 and <span style="color: red;">C000-1216</span>	C000-1115 <span style="color: red;">C000-1217</span>
	<b>IV Pole</b> Stainless steel with four bottle hooks Range: 1070 -1660 mm Tube diameter: 25 / 18 mm  <b>IV Pole Clamp C000-1325 is required.</b>	C000-1324





## Aspect 100UCPlus Optional and Replacement Accessories:

Photo	Description <b>Red text = Metric accessory rails</b>	Catalog Number
	<b>IV Pole Clamp</b> Aluminum. 25 mm pole opening	C000-1325
	<b>Clamp-on Accessory Rail</b> Rail is 8 1/2" long x 3/8" thick x 1-1/8" high. <b>Clamp-on Accessory Rail</b> 10mm x 25mm x 216mm clamp-on accessory rail.	A100-2218 A100-2838
	<b>Stirrups Mounting Rail extension Kit</b> Included with table. For use with leg-holders (Leg-holders are sold separately).	A100-1876 A100-2676
	<b>Hand Control for 100UCPlus and 100UC.</b> Controls height up/down, trendelenburg/reverse, slide, cross table movements. (Included with table).	K000-0158
	<b>Foot Switch Control for 100UCPlus.</b> (Included with table).	A100-1540
	<b>Side-rail clamp</b> Rotating, with socket for attachment of accessories with 5/8 inch (15.88mm) circular post. Requires A100-1007 (or A100-2379 for metric rails) clamp-on accessory rail.	C000-0746 C000-1065
	<b>Removable radiographic cassette tray.</b> Loads from either side of table.	A100-3139

## Aspect 100UCPlus Consumables

Photo	Description <b>Red text = Metric accessory rails</b>	Catalog Number
	<b>Vinyl Urology Table Covers</b> Box of 20	C000-0645
	<b>Uro Drain Bags with Hose</b> Disposable drain bag and hose assembly. Box of 20	<b>Sterile:</b> # C000-0593  <b>Non-sterile:</b> #C000-1111
	<b>Disposable Covers for Aspect 100UCPlus Foot Controls</b> 24" x 20" Box of 50	C000-0492

## Aspect 100UCPlus Consumables

Photo	Description <b>Red text = Metric accessory rails</b>	Catalog Number
	<b>Uro Drain Support Hoop</b> Collapsible & removable drain bag support hoop.	C000-0595
	<b>3\" data-bbox="351 235 718 311"/&gt;</b>	C000-0571
	<b>Uro Collection Unit</b> 5 gallon/18.9 liter capacity capacity, Box of 10. Soft sided Vinyl container gently expands into a cube as it collects waste fluid when floor drain or suction is unavailable. Internal baffle system enables container to fill evenly as powder converts liquid to gel. Built in zipper and handles allow for storage, closure of drain bag and hose assembly for enhanced infection control during transport and disposal.	C000-0612
	<b>Uro Collection Unit</b> 15 liter capacity, Box of 10. Flexible plastic sided container expands into a cube as it collects waste fluid when floor drain or suction is unavailable. Internal baffle system enables container to fill evenly as powder converts liquid to gel. Integrated valve system allows user to shut off the flow of liquid from the drain bag in order to change collection units..	C000-1307

## Aspect 100UCPlus Shipping Crate Specifications

- Weight (approx.): 815 lbs. (370kg)
- Length: 96" (2438mm)
- Width: 42" (1067mm)
- Height: 47" (1194mm)



## Aspect 100-US Gyn/Fertility Clinic Table



**Aspect 100US A100-1479 (120VAC)    A100-2674 (230VAC)**

GYN/Fertility Mobile Imaging Table.




Table elevation: 26" - 38" (66 - 96.5cm)

An economical C-arm table for hysterosalpingograms.


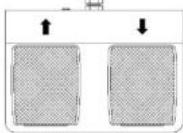



Includes the following:

- A100-758 Foot control
- X100-1687: tabletop pad, 24" x 46" (610mm x 1168mm) 1" (25.4mm) thick
- C000-0571: Patient restraint straps
- A100-2380: Leg support tabletop extension with pad
- A100-2248: Basic Knee Crutches [pair] with a pair of rotating side rail clamps [C000-0746]

### Aspect 100US Optional and Replacement Accessories:

Photo	Description <b>Red text = Metric accessory rails</b>	Catalog Number
	<b>Replacement Tabletop Pad</b> 2 inch thick foam pad. 24" wide x 46" long.	X100-1687
	<b>Armboard</b> Quick-release, rail mounted.	A310-056 <b>A310-062</b>
	<b>Replacement pad</b> for Armboard A310-056 or A310-062	X100-1755
	<b>Hand control for 100US</b>	K000-0041

## Aspect 100US Optional and Replacement Accessories:

Photo	Description <b>Red text = Metric accessory rails</b>	Catalog Number
	<b>Basic Knee Crutches</b> For patients up to 400 lbs. Basic pair with pads. Includes two C000-0746 rotating side rail clamps.	A100-2248 <b>A100-2677</b>
	<b>Foot Control 100US</b> (Included with table.)	A100-758
	<b>Side-rail clamp</b> Rotating, with socket for attachment of accessories with 5/8 inch (15.88mm) circular post. Requires A100-1007 (or A100-2379 for metric rails) clamp-on accessory rail.	C000-0746 <b>C000-1065</b>
	<b>IV Pole</b> Stainless steel with four bottle hooks Range: 1070 -1660 mm Tube diameter: 25 / 18 mm  <b>IV Pole Clamp C000-1325 is required.</b>	C000-1324
	<b>IV Pole Clamp</b> Aluminum. 25 mm pole opening.	C000-1325

## Aspect 100US Shipping Crate Specifications

Weight (approx.): 245 kg / 540 lbs.  
Length: 244 cm / 96 in. Width: 107 cm.  
42 in. Height: 119 cm / 47 in.

## MDS Video Integration System



### MDS Description # A280-0330 (120VAC) #A280-0335 (230 Volt)

MDS Mobile Display System, Complete with High Definition Video Monitors, Router and tablet controller.

- 120VAC is standard.
- 230VAC is optional and must be specified at time of order

Includes the following:

- Dual-Arm Mobile Display Stand with adjustable counterbalanced arms
- Adjustable shelves for placement of ancillary equipment (Qty. of 4)
- Bilateral input source cable manifold supporting DVI, SDI, VGA, S-Video, SOG and composite inputs
- Built-in isolation transformer and 8 outlet power strip for video and ancillary equipment AC power
- FSN 26" (660mm) HD LED Widescreen Multi-Modality Video Monitors (Qty of 2).
- FSN IPS1000 Video Router/Controller with 10.1" Tablet Controller and articulated support arm with lock
- Dual 8" (20cm) "Preview" LCD Reference Monitors, mounted to MDS stand center column
- Qty (1) External cable storage hook
- Qty (1) Kit - Accessory Basket Kit
- Protective covers with storage for 26" (660mm) FSN monitors (Qty of 2)
- Video Cable kit, including common cables/connectors needed for connection of C-arm and equipment that will be mounted onto the MDS Mobile Display System. Customer choice of one of the following cable kits is included with the MDS System.





**Customers must specify cable kit choice with order:**

- **A280-0325:** Cable kit for GE OEC 9800/9900 C-Arms and camera cable kit for systems mounted onto the MDS
- **A280-0328:** Cable kit for all DVI output C-Arms and camera cable kit for systems mounted onto the MDS
- **A280-0326:** Cable kit for Siemens VGA C-Arm with camera cable kit for systems mounted onto the MDS
- **A280-0425:** Cable kit for Siemens Cios C-Arm with Video Manager and camera cable kit for systems mounted on the MDS

## MDS Optional Cable Kits and Replacement Accessories:

Photo	Description	Catalog Number
	<b>External Source Cable Kit</b> Add this Kit if the customer plans to add any external device that is not mounted directly onto the MDS stand other than the C Arm (i.e. ultrasound, PACS, hemodynamics) Kit Includes: qty (1) each of 20 ft. (6.1m) S Video, (1) 23 ft. (7m) HDMI, 20 ft. (6.1m) SVG cable and qty. (2) short DVI to HDMI adapters.	A280-0437
	<b>Endo Shelf</b> Bilateral slide-out shelf, increases storage and shelf space For MDS or MDS D.	A280-0415
	<b>CO<sub>2</sub> gas brackets</b> Pair of holders for type E gas cylinders (for insufflator use) mounted to sides of MDS stand.	A280-0401
	<b>Medical Rail Shelf</b> Shelf Includes medical rails. For accessory attachments.	A280-0130
	<b>Shelf with Keyboard</b> Kit includes sliding keyboard carriage. Suitable for most standard customer-supplied keyboards.  Will attach to existing shelves.	A280-0526
	<b>Accessory Basket Kit</b> For customer-supplied foot switch or other storage.	A282-0007
	<b>Suction Canister Kit</b> Includes: Suction canister rail clamps, bracket and 3000 ml canister.	A282-0009
	<b>Cable Hook</b> Cable storage J-hook for external cables.	A280-0570

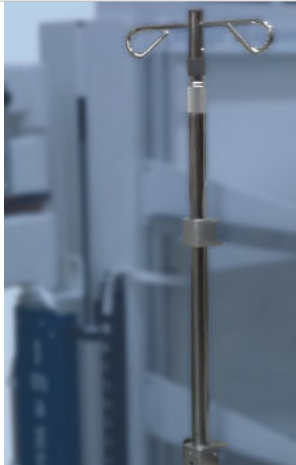
## MDS Optional and Replacement Accessories:

Photo	Description	Catalog Number
	<b>IV Pole Kit</b> Includes two IV rail clamps suitable for mounting a 4' 2-hook IV pole onto two medical rail equipped (A280-0130) shelves. Note: this option requires two shelves with medical rails. This IV pole kit is for MDS systems only.	A280-0508
	<b>MDS Monitor Protection Covers</b> For 26" (660mm) flat panel monitor. (2 covers are included with MDS System A280-0330)	X280-0262
	<b>Replacement IPS1000A Control OR</b> Video processing system used in conjunction with MDS to scale, reformat and split video signals into both monitors.	<b>C000-1041</b> GE OEC Siemens Arcadia  <b>C000-1140</b> GE Elite Siemens Cios Ziehm RFD Philips
	<b>Replacement ITouch Screen Tablet with dock</b>	C000-1043

## MDS Installation Options

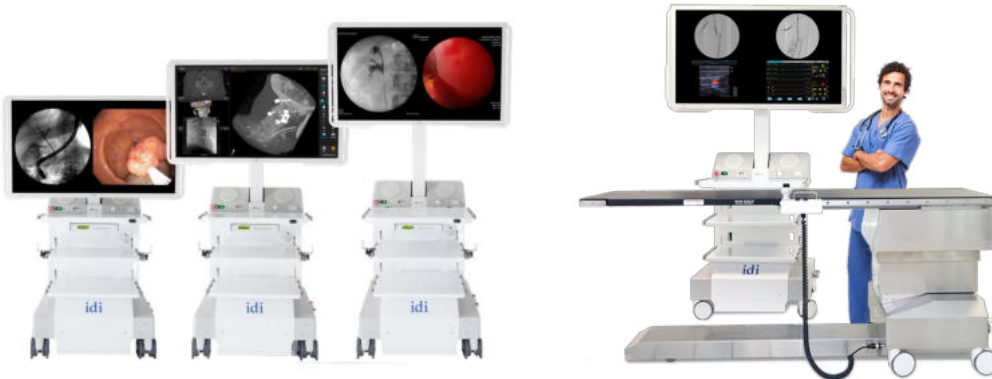Description	Catalog Number
<b>MDS Installation options</b> System installation by IDI, including two-day operator in-service training, is available for the USA and Canada. Dealers who choose to have IDI perform the technical installation only will receive a discount provided they supply two days of certified clinical training. Additional in-service training of more than two days is available on a billable basis only. Contact IDI regarding fees for additional training. Installation costs are not discountable and installation is only available in the USA and Canada. Installation outside of the USA and Canada is the responsibility of the dealer.	
<b>MDS Installation, technical integration and clinical training by IDI in the continental United States</b>	Y-MDSINS-01
<b>MDS installation and integration only by IDI technician in Continental United States</b> OEM agrees to supply 2-day certified clinical applications training to customer.	Y-MDSINS-04
<b>MDS Installation, integration and clinical training by IDI in Alaska or Hawaii</b> 4-day installation, integration and clinical operation training by IDI.	Y-MDSINS-02
<b>MDS Installation, integration and clinical training by IDI in Canada</b> 4-day installation, integration and clinical operation training by IDI.	Y-MDSINS-03

## ilex55 Video Integration System



### ilex55 Description

The Ilex 55 mobile video platform features a 55"/139cm 4K UHD resolution monitor on a mobile platform. The system can display one 4K UHD or multiple HD resolution clinical images in a single, split screen, or quad screen mode.

#### ilex55

#A280-0600 (120VAC)

#A280-0781 (230 Volt)

Includes the following:

- 55"/139cm 4K UHD medical monitor
- Three equipment shelves 19" (48cm) x 28" (71cm) and one platform 19" (48cm) x 20" (51cm)
- Two video cable storage hooks
- 5"/127mm casters with lock, free float and directional control
- Support column with powered 20" (508mm) vertical travel and 0°, 2.5° and 5° manual tilt
- Bluetooth speaker system with volume control
- Configurable video signal input panels with the following available inputs: 4x DVI-D, 2x DisplayPort, 2x VGA, 1x S-Video, 1x SDI & 2x SOG.
- Four USB ports
- Two 120v power ports
- Remote control capability for monitor PIP functions: single, side-by-side or quad modes
- Vertical transport handles
- Cable kit suitable for connecting one SOG Lex/Right C-Arm, one DVI Left/Right C-Arm, one VGA Ultrasound or hemodynamic monitor, and one DVI ultrasound or camera.
- Two boxes of sterile remote control drapes 15 per box

## ilex55 Power Options

### Option 1: North American Version, 15 amp.

For connection to 110-120 VAC, 60Hz, 15 amp outlet, 12 amp max load rating. Includes one built-in 800 VA isolation transformer and 20 ft., 20 amp power cord, 8 position power strip (12A max draw)



## ilex55 Dimensions and weight

- Base cover: 31in.(79cm) wide x 29in.(74cm) deep.
- Maximum base footprint: 36in.(92cm) wide x 35in.(89cm) deep.
- Monitor Handles are 55in.(140cm) wide.
- Height with monitor retracted is 72in.(182cm)
- Height with monitor extended is 92in.(234cm)
- Weight: with monitor 460 lbs. (209 kg), without 360lbs (163 kg)

## ilex55 Monitor

- Screen Size: 55" medical grade TFT AM LCD/IPS Technology/LED backlight
- Resolution: 4K UHD resolution 3840x2160 pixels 16:9 Aspect
- Brightness: Max: 500 cd/m<sup>2</sup> (typ)
- Response: Gray-to-gray = 8 msec (typ)
- Weight: 53 lbs. (24kg)
- Dimensions: 49.6" (126cm) W x 28.9" (73.4cm) H x 3.4" (8.7cm) D

## ilex 55 optional and Replacement Accessories

	<b>Protective Monitor and Cart Covers</b> Padded cloth cover protects the monitor. Matching skirt covers the cart's shelves and main frame. Skirt comes with buckles and straps in the back.	<b>C280-0805</b> Set <b>C280-0817</b> Monitor cover <b>C280-0866</b> Cart cover
	<b>23" Monitor with Magnetic Stand</b> Additional 23" monitor for use with ilex55 mounted on a magnetic base and tillable for maximum visibility.	A280-0977

## ilex55 Installation Options

Installation by IDI in continental US. Includes one-day operator training.	Y-ILEXINS01
Installation by IDI in Alaska or Hawaii, including one-day operator training.	Y-ILEXINS02



## ilex32 4K video monitor and stand



### ilex32 4K Description

The Ilex32 4K monitor and stand provides a cost effective solution for those customers requiring an additional UHD 4K monitor or monitors with mirrored images in multiple locations. Designed for use as a large table side monitor, the ilex 32 4K combines practical features with easy use.

### ilex32 4K # A280-0821 120/230 VAC)

Includes the following:

- 31.5" UHD 4K display features signal input options for HDMI, DVI, DisplayPort plus DC out power for smaller components
- Single or split-screen images
- Heavy duty cart stand with tilt (Maximum height: 76"/193cm to top of monitor)
- Cable channel
- Accessory basket
- 15 foot power cord (For use outside the USA please specify plug type)
- Cable hook
- Qty. (1) C-Arm kit including a pair of input cables qty 2
- One year parts and labor warranty

### ilex32 4K Installation options

The customer typically performs installation. Installation is not available outside of the United States. For Dealers who choose to have IDI perform the installation:

#### ilex32 4K Installation by IDI in the continental US

Installation by IDI in continental US. Includes one-day operator training.

Y-ILEXINSO1

#### ilex32 4K Installation by IDI in Alaska or Hawaii

Installation by IDI in Alaska or Hawaii, including one-day operator training.

Y-ILEXINSO2

## ilex32 video monitor and stand



### ilex32 Description

The ilex32 monitor and stand provides a cost effective solution for those customers requiring an additional monitor or monitors with mirrored images in multiple locations.

- 120VAC is standard
- 230VAC is optional and must be specified at time of order

### ilex32 # A280-0630 (120VAC)    #A280-0743 (230 Volt)

Includes the following:

- 32"/81cm HD medical grade monitor with medical bios and input options for Dual SOG, Dual DVI, VGA, SDI, S Video
- Single or split-screen images
- Heavy duty cart stand with tilt (Maximum height: 76"/193cm to top of monitor)
- Cable channel
- Accessory basket
- 15 foot power cord
- Cable hook
- Qty. (1) C-Arm kit including a pair of input cables qty 2
- One year parts and labor warranty

### ilex32 Installation options

The customer typically performs installation. Installation is not available outside of the United States. For Dealers who choose to have IDI perform the installation:

<b>ilex32 Installation by IDI in the continental US</b> Installation by IDI in continental US. Includes one-day operator training..	Y-ILEXINSO1
<b>ilex32 Installation by IDI in Alaska or Hawaii</b> Installation by IDI in Alaska or Hawaii, including one-day operator training.	Y-ILEXINSO2

## ilex C Utility Cart With Accessory Rail



Shown with motorized  
C-Arm control panels

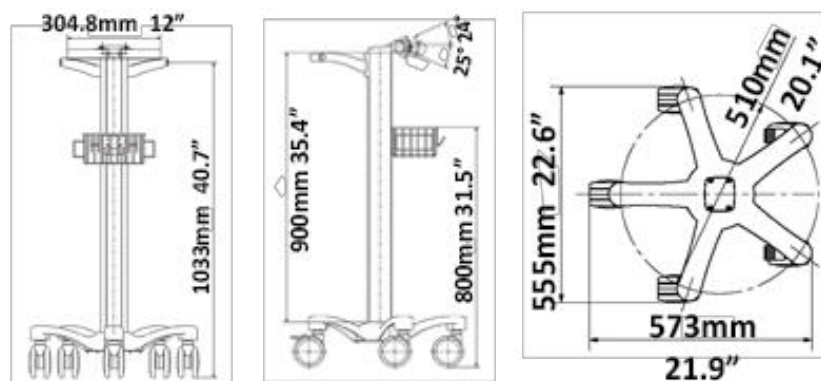
### ilex C Description # A280-0893

The IlexC is a utility cart for mounting motorized C-Arm control panels, patient monitoring and other accessory rail-mounted devices. Load capacity: 30 lbs. (13.6kg) Complete with adjustable U.S. and metric 12" (305MM) medical accessory rail and accessory basket.

- Two locking casters
- Utility basket with cable wraps
- 12" (305mm) adjustable medical accessory rail (U.S.  $\frac{3}{8}$ " x  $1\frac{1}{8}$ " & metric 10mm x 25mm)
- Accessory rail is 42" (1067mm) from the floor
- Control handle
- Up to 30 lbs (13 kg) capacity
- Approximately 22" x 23" footprint

### ilex C Dimensions

Note: Mounting rail is metric: 10mm x 25mm



## Connexus Direct Connect video transmission

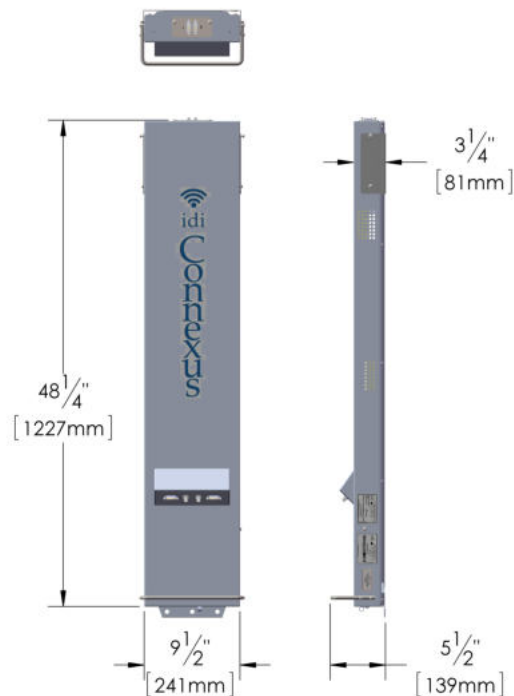


### Connexus Direct Connect Description # A280-0916

Wall mounted signal conversion system for up to two channels of uncompressed, encrypted HD video. Attach C-Arm control cart video cables to the Connexus Direct Connect via DVI-I connectors on the top of the unit to transmit left and right channel C-Arm signals to any location in a room; wall monitor, boom monitors or IDI's Ilex55, Ilex32 4K or Ilex 32 video integration systems.

### Connexus Direct Connect Dimensions/Weight

Weight: 17 lbs (7.71kg)



## SIMEON Lighting Systems

### Sim.LED 4500 Mobile Light



### Sim.LED 4500 Mobile Light Description #K400-0103

The **Sim.LED 4500** is the ideal light for outpatient and plastic surgery. This light provides ideal illumination, uniform light and very good depth illumination.

### Sim.LED 4500 Mobile Light Technical Data

Central illumination intensity at a distance of 1 meter [lx]	100,000
Electrical dimmability from/to [lx]	25,000 – 100,000
Color temperature [K]	4,500
Field size d10 at a distance of 1 meter [mm]	170
Color rendering index Ra	96
Red rendering index R9	96
Illumination depth without refocusing (L1/L2) [mm]	1,530
Light head weight [kg]	3
Light head power consumption [W]	28
No. of LEDs	18
Approval	CE
LED lifespan [h]	> 60,000

Tolerance  $\pm 10\%$

Technical specifications are subject to change

## More SIMEON Lighting Systems

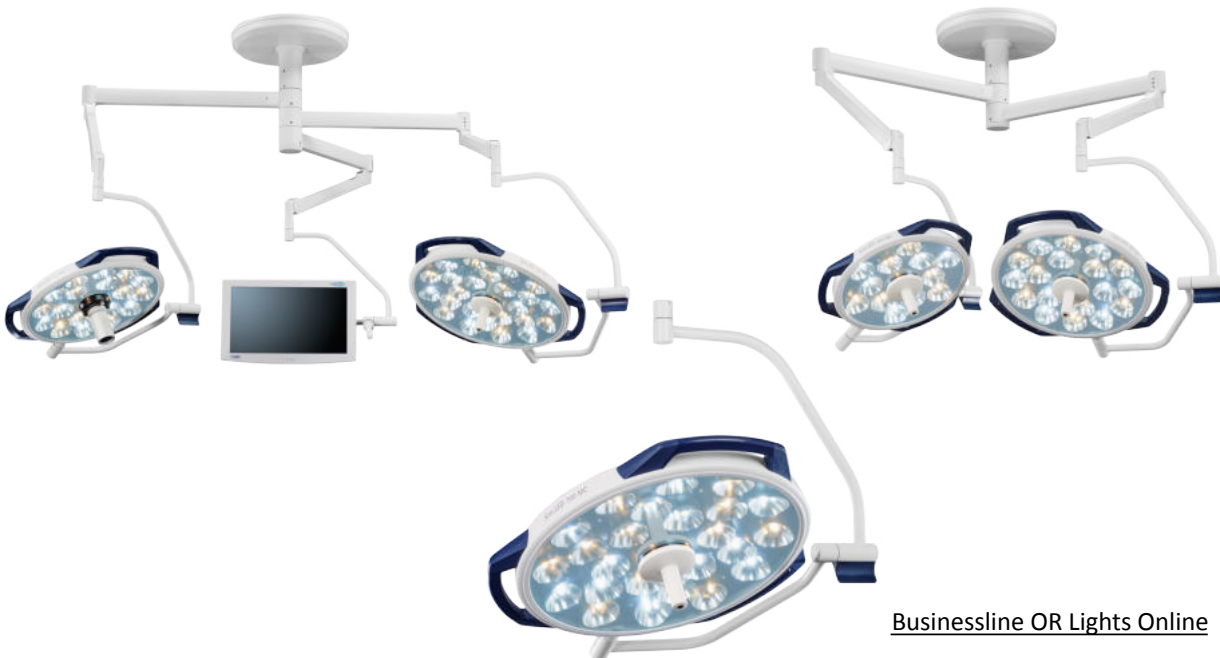
IDI offers a wide variety of premium O.R. Light options from SIMEON Medical.  
[CLICK HERE](#) to see all of IDI's SIMEON lighting systems.

### SIMEON Highline Lighting Systems Examples



[Highline OR Lights Online](#)

### SIMEON Businessline Lighting Systems Examples



[Businessline OR Lights Online](#)

## Aspect Four-Way Float Table Specifications & Accessories

•=Standard feature ◊=Optional NA=Not available	100 - 4	100 - 4T G3	ISR G3
Elevate (Height with 2" pad)	29.9 - 39.7" 76 - 101cm	32-45" 81-114 cm	33-45" 84-114 cm
Trendelenburg Tilt	NA	±12°	±12°
Lateral Tilt	NA	NA	±12°
Longitudinal Tabletop Motion	40" float 102 cm	32" float 81 cm	32" float 81 cm
Transverse Tabletop Motion	10" float 25 cm	8" float 20 cm	8" float 20 cm
Standard Tabletop Size	21/30" x 80" 53/76 x 203 cm	24"/22"/19.5" x 84" 61/56/50 x 213cm	24"/22"/19.5" x 84" 61/56/50 x 213cm
Imaging Area	72" 183 cm	68" / 173cm  80" / 203cm with optional peripheral extension	68" / 173cm  80" / 203cm with optional peripheral extension
Patient Capacity	450 lb./204 kg	600 lb./272 kg	600 lb./272 kg
Table Weight	490 lb. 222 kg	800 lb. 363 kg	800 lb. 363 kg
Tabletop Pad	•	•	•
Restraint Straps	•	•	•
Pendant Hand Control	NA	•	•
4-way Panning Hand Control	•	•	•
Foot Control	NA	◊	◊
Battery operation	Backup for elevation only	• Up to 5 days Will vary based on case load and patient size	• Up to 5 days Will vary based on case load and patient size
Armboard, rail mounted	NA	◊	◊
Armboard, carbon fiber	◊	◊	◊
VAB Armboard	◊	◊	◊
Clamp-on Accessory Rail	◊	◊	◊
Catheter Tray extension	◊	◊	◊
Peripheral Extension	NA	◊	◊
Headrest extension	NA	◊	◊
Table-Mounted Radiation Shield	◊	◊	◊
Anesthesia Screen Holder	NA	◊	◊
IV Pole	◊	◊	◊

## Aspect Urology & GYN Table Specifications & Accessories

<b>•Standard feature</b> <b>◊Optional NA Not available 100 UC-Plus</b>	<b>100 UC-Plus</b>	<b>100 US</b>
<b>Elevate (Height with pad)</b>	33.3" – 43.1" (84.5 – 109.4 cm)	27" – 39" (68.5 – 99cm)
<b>Trendelenburg Tilt</b>	± 15°	NA
<b>Transverse travel</b>	± 3.5", 7" total (± 8.9 cm, 17.8 cm total)	NA
<b>Longitudinal Tabletop Motion</b>	10" (25.4 cm)	NA
<b>Tabletop Size</b>	28"W x 46"L (71 x 117 cm)	24"W x 46"L (61 x 117 cm)
<b>Tabletop Extension Size</b>	24"W x 34"L (61 x 86 cm)	NA
<b>Imaging Area</b>	Main tabletop: 25"W x 20"L ( 63.5 x 50.8 cm)	22"W x 32"L (56 x 81 cm)
<b>Imaging Area with table extension</b>	54" (137 cm)	NA
<b>Patient Weight Capacity</b>	Upper body on main table: 525lb. (238 kg) Upper body on extension: 400 lbs ( 181kg) Load capacity on extension only: 200 lbs (90.7kg)	450 lbs (204 kg)
<b>Table Weight</b>	490 lb. (222 kg)	375 lbs ((170kg)
<b>Tabletop Pad (2" thick)</b>	•	•
<b>Tabletop extension with pad</b>	•	Leg support tabletop extension
<b>Hand control pendant &amp; footswitch control</b>	•	Footswitch only
<b>Patient restraint straps</b>	•	•
<b>Removeable drain bag support hoop</b>	•	NA
<b>Disposable drain bag &amp; hose assembly, box of 20</b>	•	NA
<b>Removable radiation shield flap for perineal end of table</b>	•	NA
<b>Armboard, rail-mounted with pad and strap</b>	•	◊
<b>Split Leg Stirrups</b>	◊	NA
<b>Disposable covers for footswitch, box of 50</b>	•	◊
<b>Radiographic cassette tray</b>	◊	NA



## Aspect Pain Management Table Specifications & Accessories

•Standard feature ◊Optional NA Not available	100E	100T	100RT	100RTL
Elevate (Height with 2" pad)	28 - 39.8" 71 - 101 cm	30 - 41.8" 76 - 106 cm	32.4 - 44.2" 82 - 112 cm	33.2 - 43.2" 84.3 - 109.7 cm
Trendelenburg Tilt	NA	±15°	±15°	±15°
Lateral Tilt	NA	NA	20°	20°
Longitudinal Tabletop Motion	NA	NA	NA	10"
Standard Tabletop Size	24" x 80" 61 x 203 cm	24" x 80" 61 x 203 cm	24" x 80" 61 x 203 cm	24" x 80" 61 x 203 cm
Optional Tabletop Size	22" x 80" 56 x 203 cm with facial cut-out	22" x 80" 56 x 203 cm with facial cut-out	22" x 80" 56 x 203 cm with facial cut-out	22" x 80" 56 x 203 cm with facial cut-out
Optional GI Tabletop Size	28" x 80" 71 x 203 cm	28" x 80" 71 x 203 cm	28" x 80" 71 x 203 cm	28" X 80" 71 x 203 cm
Imaging Area	66" 167.6 cm	62" 157.5 cm	62" 157.5 cm	59" 149.9 cm
Patient Capacity	450 lb. 204 kg	450 lb. 204 kg	450 lb. 204 kg	450 lbs. 204 kg
Table Weight	290 lb. 131 kg	325 lb. 147 kg	380 lb. 172 kg	355 lbs. 161 kg
Tabletop Pad (2" thick)	•	•	•	•
Restraint Straps	•	•	•	•
Facial Support Cushion	•	•	•	•
Pendant Hand Control	•	•	•	•
Foot Control	◊	◊	◊	◊
Armboard, rail mounted	◊	◊	◊	◊
Armboard, carbon fiber	◊	◊	◊	◊
VAB Armboard	◊	◊	◊	◊
Clamp-on Accessory Rail	◊	◊	◊	◊
Scatter-Guard Radiation Shield	◊	◊	◊	◊

## Image Diagnostics, Inc.

Since 1986, IDI has been a manufacturer of innovative products for surgical, interventional and diagnostic imaging.

Our product focus is in electro-mechanical solutions for fluoroscopic and digital radiographic applications.

IDI is committed to superior product design, value and customer support in everything we do.

### Image Diagnostics, Inc.

310 Authority Drive  
Fitchburg, Massachusetts  
01420

For customer service:  
[sales@imagediagnostics.com](mailto:sales@imagediagnostics.com)

Toll free: 877.304.5434  
Main: 978.829.0009  
Fax: 978.829.0027

[ImageDiagnostics.com](http://ImageDiagnostics.com)

**idi**  
Image Diagnostics Inc.



A Veteran-Owned Business