

Urology and GYN/Fertility Tables

Aspect

# 100-UCPLUS 100US

Tables & Accessories



100UCPLUS Urology Table



100US GYN / Fertility Table

INTELLIGENT DESIGN ► TANGIBLE OUTCOMES

**idi**  
imagediagnostics.com

## Aspect 100-UCPLUS Urology/Multi-purpose Table



Shown with optional leg holders (#A100-2247).



Shown with tabletop extension (A100-1400)

### Aspect 100UCPlus description:

Imaging table for urology and multi-purpose use. Urological accessories are easily removable so the table can be quickly configured for multi-purpose C-arm imaging or general patient positioning for non-urological procedures.

- 120VAC is standard
- 230VAC is optional and must be specified at time of order

### Aspect 100UCPlus #A100-1517 (120VAC with U.S. Accessory Rails) A100-2639 (230VAC with metric accessory rails)

Includes the following:

- A100-1876: Stirrups Mounting Rail extension Kit
- X100-1818: tabletop pad, 28" (711mm) x 46 1/2" (1181mm) x 2" (51mm) thick
- A100-1400: Tabletop extension, 24" (610mm)W x 34" (864mm) (Metric tables A100-2675)
- X100-1817: Tabletop extension pad, 24" (610mm) x 34" (864mm) x 2" (51mm) thick
- K000-0158: Pendant hand control
- A100-1540: Foot switch control
- U100-3025: Removable x-ray shield flap, for perineal end of table
- C000-0571: Patient restraint straps
- A310-056: Arm board, rail mounted, with pad and strap (A310-062 for metric rails)
- C000-0595: Urology drain bag support hoop
- C000-0593: Disposable drain bag & hose assembly, sterile, box of 20.
- C000-0492: Disposable covers for foot switch, box of 50
- C000-0645: Disposable tabletop pad protector (qty. of 1)
- C000-0612: 5 gallon collection bag (qty of 1)

## **Aspect 100UCPlus Package Special #A100-2521 (120VAC)**

**A100-3180 (230VAC)**

Includes the following:

- A100-1517: Aspect 100UCPlus Urology Table with standard accessories
- A100-2246: Two ComfortFin® 350 stirrups with clamps for US rails
- A310-056: One rail-mounted arm board
- C000-0593: Disposable drain bag & hose assembly, sterile, box of 20
- C100-3044: One ComfortFin stirrup storage cart

## **Aspect 100UCPlus Package Special #A100-3107 (120VAC)**




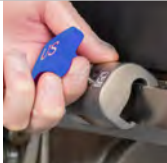

**A100-3181 (230VAC)**

Includes the following:








- A100-1517: Aspect 100UCPlus Urology Table with standard accessories
- A100-2247: Two ComfortFin 500 stirrups with clamps for US rails
- A310-056: One rail-mounted arm board
- C000-0593: Disposable drain bag & hose assembly, sterile, box of 20
- C100-3044: One ComfortFin stirrup storage cart

ComfortFin® is a registered trademark of Image Diagnostics, Inc.

## Aspect 100UCPlus Optional and Replacement Accessories:

Photo	Description <span style="color: red;">Red text = Metric accessory rails</span>	Catalog Number
	<p><b>IDI ComfortFin 350 Lithotomy Leg holders (Stirrups pair only)</b>                      ComfortFin Lithotomy Leg holders are designed to improve infection control with durable, easy to clean silicone straps. Includes integrated lateral fin to help protect the head of the fibula and the peroneal nerve. 4" more working space with an increased range of motion. Patient weight capacity: 350 lbs (159 kg).</p>	C100-3042
	<p><b>Replacement Boot with Straps</b> for IDI C100-3042 ComfortFin 350 Lithotomy Leg Holders</p>	C100-3080
	<p><b>Kit:</b> IDI ComfortFin350 Lithotomy Leg holders C100-3042 with C100-3046 economy side-rail clamps</p>	A100-2246 <span style="color: red;">A100-3074</span>
	<p><b>IDI ComfortFin 500 Lithotomy Leg holders (Stirrups pair only)</b>                      ComfortFin Lithotomy Leg holders are designed to improve infection control with durable, easy to clean silicone straps. Includes integrated lateral fin to help protect the head of the fibula and the peroneal nerve. Boots are sized for today's heavier patients. 4" (102mm) more working space with an increased range of motion. Patient weight capacity: 500 lbs. (227kg).</p>	C100-3045
	<p><b>Replacement Boot with Straps</b> for IDI C100-3045 ComfortFin 500 Lithotomy Leg Holders</p>	C100-3079
	<p><b>Kit:</b> IDI ComfortFin 500 Lithotomy Leg holders C100-3045 with C100-3046 economy side-rail clamps.</p>	A100-2247 <span style="color: red;">A100-3075</span>
	<p><b>Side-rail clamp, removable</b>                      Attaches and locks anywhere on the side rail for fast and easy mounting . A push button safely releases rail clamp —helps to prevent false rail engagement. Autoclavable.</p>	C100-3043  <span style="color: red;">C100-3076</span>
	<p><b>Economy side-rail clamp.</b> Attaches anywhere on the table rail for fast and easy mounting. Rounded ergonomic shape is easy on the hands. Autoclavable.</p>	C100-3046  <span style="color: red;">C100-3077</span>
	<p><b>Lithotomy Leg Holder Cart</b>                      Easy to push/pull with ergonomic handles. Stirrups are isolated to prevent damage. Clamp storage compartment. Designed to safely store ComfortFin stirrups when not in use.</p>	C100-3044

## Aspect 100UCPlus Optional and Replacement Accessories:

Photo	Description <span style="color: red;">Red text = Metric accessory rails</span>	Catalog Number
	<b>Replacement Tabletop Pad</b> 2" (51mm) thick foam pad. 28" (711mm) wide x 46 1/2" (1181mm) long.	X100-1818
	<b>Replacement Extension Tabletop Pad</b> 2" (51mm) thick foam pad. 28" (71mm) wide x 34" (864mm) long.	X100-1817
	<b>Single Use Stirrup Boot Cover</b> Waterproof single use stirrup boot covers help eliminate strike through and bacterial migration. Infectious patients (known or suspected). Patients with open sores, abrasion or lesions. Protect against cross contamination.	C100-3070
	<b>Physician Elbow Supports (pair)</b> Removable	C000-0594
	<b>Armboard</b> Quick-release, rail mounted (requires Clamp-on Accessory Rail A100-1007 or A100-2379 for use on tabletop extension.	A310-056 <span style="color: red;">A310-062</span>
	<b>Replacement pad for A310-056 or A310-062 Armboard</b>	C100-3329
	<b>Single use soft armband straps, package of 12</b>	C000-1555
	<b>Split Leg Holders</b> For use with Aspect 100UC Plus Tables. Includes pads. Prone percutaneous, lateral or modified lateral nephrostomy. Includes two split leg holder clamps # C000-1115 or <span style="color: red;">C000-1216</span> .	C000-1146 <span style="color: red;">C000-1216</span>
	<b>Right replacement pad for Spilt Leg Holder C000-1146 or <span style="color: red;">C000-1216</span>.</b>	C000-1357
	<b>Left replacement pad for Spilt Leg Holder C000-1146 or <span style="color: red;">C000-1216</span>.</b>	C000-1358
	<b>Replacement Split Leg Holder Clamp</b> For use with split leg holders. # C000-1146 and <span style="color: red;">C000-1216</span>	C000-1115 <span style="color: red;">C000-1217</span>
	<b>IV Pole</b> Stainless steel with four bottle hooks Range: 42" (1070mm) / 65" (1660) mm Tube diameter: .98" (25mm) / .079" (18 mm) <b>IV Pole Clamp C000-1325 is included.</b>	A100-3285

## Aspect 100UCPlus Optional and Replacement Accessories:

Photo	Description <span style="color: red;">Red text = Metric accessory rails</span>	Catalog Number
	<p><b>Clamp-on Accessory Rail</b>  <b>(For use on tabletop extension)</b>                      Rail is 8 ½" long x ¾" thick x 1-1/8" high.  <span style="color: red;"><b>Clamp-on Accessory Rail</b></span>  <span style="color: red;">Metric rail is 10mm x 25mm x 216mm</span></p>	<p>A100-2218   <span style="color: red;">A100-2838</span></p>
	<p><b>Stirrups Mounting Rail extension Kit</b>                      Included with table. For use with leg-holders                      (Leg-holders are sold separately).</p>	<p>A100-1876   <span style="color: red;">A100-2676</span></p>
	<p><b>Hand Control for 100UCPlus and 100UC.</b> Controls height up/down, trendelenburg/reverse, slide, cross table movements. (Included with table).</p>	<p>K000-0158</p>
	<p><b>Foot Switch Control for 100UCPlus.</b>                      (Included with table).</p>	<p>A100-1540</p>
	<p><b>Side-rail clamp</b>                      Rotating, with socket for attachment of accessories with 5/8 inch (15.9mm) circular post. Requires A100-1007 (or A100-2379 for metric rails) clamp-on accessory rail for use on tabletop extension.</p>	<p>C000-0746   <span style="color: red;">C000-1065</span></p>
	<p><b>Compass Contour Radiation Shield</b>                      Contoured ceiling suspended acrylic barrier radiation protection with a lead equivalence of .5mm. Enhance eye, head and neck protection from x-ray scatter. Shielding area is ~30" (760mm) x 23.5" (600mm) with a 10" (254mm) radius along the bottom edge of the leaded acrylic shield. Does not include ceiling suspension arm.</p>	<p>A610-0200</p>
	<p><b>Compass Frame Radiation Shield</b> Rectangular ceiling suspended acrylic barrier radiation protection with a lead equivalence of .5mm. Enhance eye, head and neck protection from x-ray scatter. Shielding area is ~30" (760mm) x 23.5" (600mm). Does not include ceiling suspension arm.</p>	<p>A610-0215</p>
	<p><b>Compass Contour and Frame Drape</b>                      Custom shaped drape for the Compass Contour and Frame Radiation Shield. Slips over the entire shield and is held in place with an elastic band. Sterile. Case of 20.</p>	<p>C610-0256</p>

## Aspect 100UCPlus Consumables

Photo	Description <span style="color: red;">Red text = Metric accessory rails</span>	Catalog Number
	<b>Vinyl Urology Table Covers</b> Box of 20.	C000-0645
	<b>Uro Drain Bags with Hose</b> Disposable drain bag and hose assembly. Box of 20.	<b>Sterile:</b> C000-0593  <b>Non-sterile:</b> C000-1111
	<b>Disposable Covers for Aspect 100UCPlus Foot Controls</b> 24" (610mm) x 20" (508mm) Box of 50.	C000-0492
	<b>Uro Drain Support Hoop</b> Collapsible & removable drain bag support hoop.	C000-0595
	<b>3" / 76mm Restraint Straps (Pair)</b> Fabric with velcro fasteners for Aspect 100UCPlus table.	C000-0571
	<b>Uro Collection Unit</b> 5 gallon/18.9 liter capacity capacity, Box of 10. Soft sided Vinyl container gently expands into a cube as it collects waste fluid when floor drain or suction is unavailable. Internal baffle system enables container to fill evenly as powder converts liquid to gel. Built in zipper and handles allow for storage, closure of drain bag and hose assembly for enhanced infection control during transport and disposal.	C000-0612
	<b>Uro Collection Unit</b> 15 liter capacity, Box of 10. Flexible plastic sided container expands into a cube as it collects waste fluid when floor drain or suction is unavailable. Internal baffle system enables container to fill evenly as powder converts liquid to gel. Integrated valve system allows user to shut off the flow of liquid from the drain bag in order to change collection units.	C000-1307

## Aspect 100UC*Plus* Shipping Crate Specifications

- Weight (approx.): 815 lbs. (370kg)
- Length: 96" (2438mm)
- Width: 42" (1067mm)
- Height: 47" (1194mm)



## Aspect 100-US Gyn/Fertility Clinic Table



### Aspect 100US A100-1479 (120VAC with U.S. accessory rails) A100-2674 (230VAC with metric accessory rails)

GYN/Fertility Mobile Imaging Table.




An economical C-arm table for hysterosalpingograms.

- 115VAC is standard.

Includes the following:

- A100-758: Foot control
- X100-1687: Tabletop pad, 24" x 46" (610mm x 1168mm) 1" (25.4mm) thick
- C000-0571: Patient restraint straps
- A100-2380: Leg support tabletop extension with pad
- A100-2248: Basic Knee Crutches [pair] with a pair of rotating side rail clamps [C000-0746]

### Aspect 100US Optional and Replacement Accessories:

Photo	Description <span style="color: red;">Red text = Metric accessory rails</span>	Catalog Number
	<b>Replacement Tabletop Pad</b> 2" (51mm) thick foam pad. 24" (610mm) wide x 46" (1168mm) long.	X100-1687
	<b>Armboard</b> Quick-release, rail mounted.	A310-056 <span style="color: red;">A310-062</span>
	<b>Replacement pad for A310-056 or A310-062 Armboard</b>	C100-3288
	<b>Single use soft armband straps, package of 12</b>	C000-1555
	<b>Hand control for 100US</b>	K000-0041

## Aspect 100US Optional and Replacement Accessories:

Photo	Description <span style="color: red;">Red text = Metric accessory rails</span>	Catalog Number
	<b>Basic Knee Crutches</b> For patients up to 400 lbs. (180kg) Basic pair with pads. Includes two C000-0746 rotating side rail clamps.	A100-2248 <span style="color: red;">A100-2677</span>
	<b>Foot Control 100US</b> (Included with table.)	A100-758
	<b>Side-rail clamp</b> Rotating, with socket for attachment of accessories with 5/8 inch (15.9mm) circular post. Requires A100-1007 (or A100-2379 for metric rails) clamp-on accessory rail for use on tabletop extension.	C000-0746 <span style="color: red;">C000-1065</span>
	<b>IV Pole</b> Stainless steel with four bottle hooks Range: 42" (1070mm) / 65" (1660) mm Tube diameter: .98" (25mm) / .079" (18 mm) <b>IV Pole Clamp C000-1325 is included.</b>	C000-1324
	<b>Compass Contour Radiation Shield</b> Contoured ceiling suspended acrylic barrier radiation protection with a lead equivalence of .5mm. Enhance eye, head and neck protection from x-ray scatter. Shielding area is ~30" (760mm) x 23.5" (600mm) with a 10" (254mm) radius along the bottom edge of the leaded acrylic shield. Does not include ceiling suspension arm.	A610-0200
	<b>Compass Frame Radiation Shield</b> Rectangular ceiling suspended acrylic barrier radiation protection with a lead equivalence of .5mm. Enhance eye, head and neck protection from x-ray scatter. Shielding area is ~30" (760mm) x 23.5" (600mm). Does not include ceiling suspension arm.	A610-0215
	<b>Compass Contour and Frame Drape</b> Custom shaped drape for the Compass Contour and Frame Radiation Shield. Slips over the entire shield and is held in place with an elastic band. Sterile. Case of 20.	C610-0256

## Aspect 100US Shipping Crate Specifications

- Weight (approx.): 245 kg / 540 lbs
- Length: 244 cm / 96 in. Width: 107 cm
- Height: 119 cm / 47 inch

## Aspect Urology & GYN Table Specifications & Accessories

• = Standard feature    ◊ = Optional    NA =Not available

	<b>100 UC-Plus</b>	<b>100 US</b>
<b>Elevate (Height with pad)</b>	33.3" – 43.1" (845mm – 1094mm)	27" – 39" (686mm – 991mm)
<b>Trendelenburg Tilt</b>	± 15°	NA
<b>Transverse travel</b>	± 3.5", 7" total (± 89mm, 178mm total)	NA
<b>Longitudinal Tabletop Motion</b>	10" (254 mm)	NA
<b>Tabletop Size</b>	Type to enter text (710 x 1170mm)	24"W x 46"L (610mm x 1168mm)
<b>Tabletop Extension Size</b>	24"W x 34"L (610 x 860mm)	NA
<b>Imaging Area</b>	Main tabletop: 25"W x 20"L (635mm x 508 mm)	22"W x 32"L (559mm x 813mm)
<b>Imaging Area with table extension</b>	54" (1370 mm)	NA
<b>Patient Weight Capacity</b>	Upper body on main table: 525lb. (238 kg) Upper body on extension: 400 lbs ( 181kg) Load capacity on extension only: 200 lbs (90.7kg)	450 lbs (204 kg)
<b>Table Weight</b>	490 lb. (222 kg)	375 lbs ((170kg)
<b>Tabletop Pad 2" (51mm) thick</b>	•	•
<b>Tabletop extension with pad</b>	•	For leg support only
<b>Hand control pendant</b>	•	•
<b>Footswitch control</b>	•	◊
<b>Patient restraint straps</b>	•	•
<b>Removeable drain bag support hoop</b>	•	NA
<b>Disposable drain bag &amp; hose assembly, box of 20</b>	•	NA
<b>Removable radiation shield flap for perineal end of table</b>	•	NA
<b>Armboard, rail-mounted with pad and strap</b>	•	◊
<b>Split Leg Stirrups</b>	◊	NA
<b>Disposable covers for footswitch, box of 50</b>	•	◊



Since 1986, IDI has been a manufacturer of innovative products for surgical, interventional and diagnostic imaging.

Our product focus is in electro-mechanical solutions for fluoroscopic and digital radiographic applications.

IDI is committed to superior product design, value and customer support in everything we do.

**Image Diagnostics, Inc.**

310 Authority Drive  
Fitchburg, Massachusetts  
01420

For customer service:  
[sales@imagediagnostics.com](mailto:sales@imagediagnostics.com)  
Telephone: 978.829.0009

**INTELLIGENT DESIGN ► TANGIBLE OUTCOMES**



Image Diagnostics, Inc. 310 Authority Drive Fitchburg, MA, 01420 USA  
A veteran-owned company since 1986. Designed and manufactured in the U.S.  
Tel: 978-829-0009 Image Diagnostics, Inc. reserves the right to update product specifications without notification. \* Warranty information for IDI products is available upon request from Image Diagnostics, Inc.

