

ASPECT 360 TABLE PACKAGE

3D Imaging without compromise

The most advanced 3D imaging table



Table available in the US through the partnership of IDI and Medifa GmbH



idi
imagediagnostics.com

INTELLIGENT DESIGN ► TANGIBLE OUTCOMES

The Aspect 360 3D Package

The first mobile, fully programmable, motorized table with joystick-control and user definable position presets.



- Unrestricted 360° imaging
- Overcome heavy patient limitations of manual tables with a fully motorized table top and fingertip panning control
- Reduce arm and hand fatigue when positioning patients during imaging - up to 551 lbs./250kg
- Joystick-based motion control of longitudinal and transverse positions
- LCD handset provides advanced motion control with three user programmable position presets
- Carbon fiber accessory rails minimize artifacts and facilitate the attachment of accessories and head support for 3D imaging of the cervical spine
- Short - Low profile base provides excellent C-arm access and positioning
- Optional width extension accessory provides an additional 4"/100mm in table width (2"/50mm per side) for positioning large patients



Aspect 360 Mobile Package Includes:

- Medifa 8000 table with motorized top (120V-230V)
- Metric accessory rails
- Two patient restraint straps
- One clamp for round rod accessories and one rail adapter

Selected optional accessories:



Width extension
24"/600mm x 2"/50mm



Radiolucent skull clamp



Joystick



Radiolucent head rest system for adults



Multi-position arm support



Footswitch



Arm shield



Armboard



Siderail clamp

Range of table top motions and dimensions

- Width = 21.25" / 540mm wide
- Length = 91.9" / 2335mm long
- 25° Trendelenburg / Reverse
- 20° roll
- Longitudinal travel 11.8" / 300mm
- Transverse travel 7" / 180mm






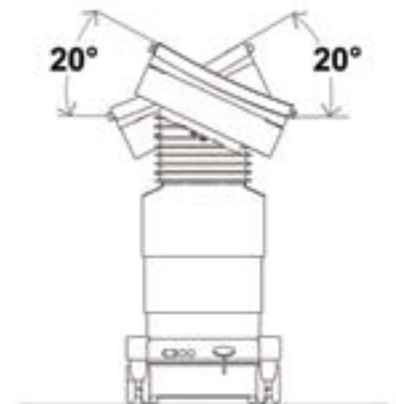


Image Diagnostics, Inc. 310 Authority Drive Fitchburg, MA 01420 Tel: 978-829-0009
 A veteran-owned company since 1986. Designed and manufactured in the U.S.
 Image Diagnostics, Inc. reserves the right to update product specifications without notification.
 Warranty information for IDI products is available upon request from Image Diagnostics, Inc.
 L100-3356 Rev. A

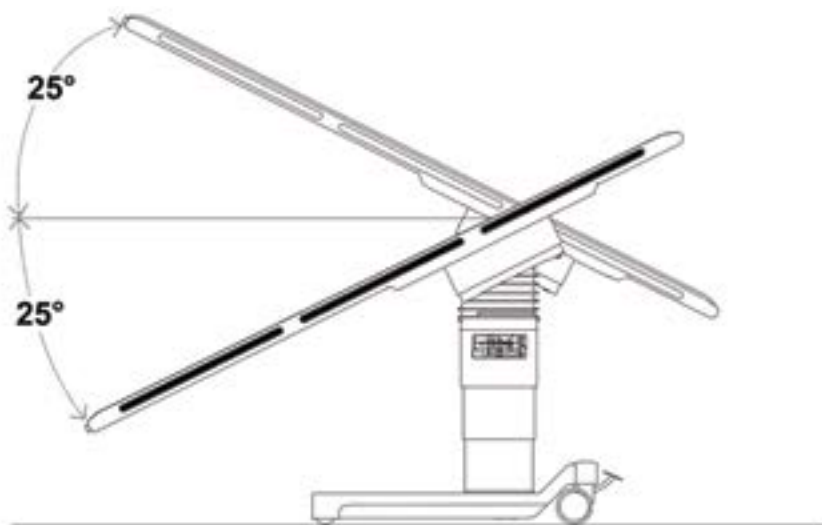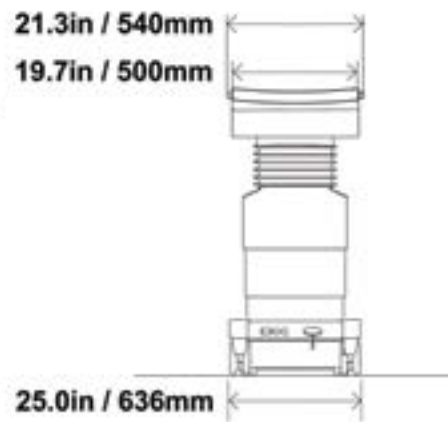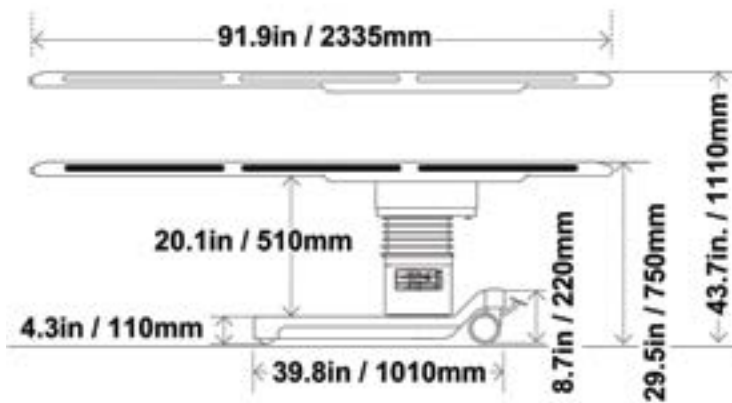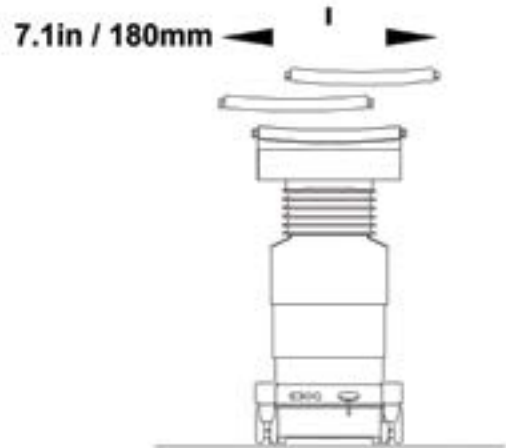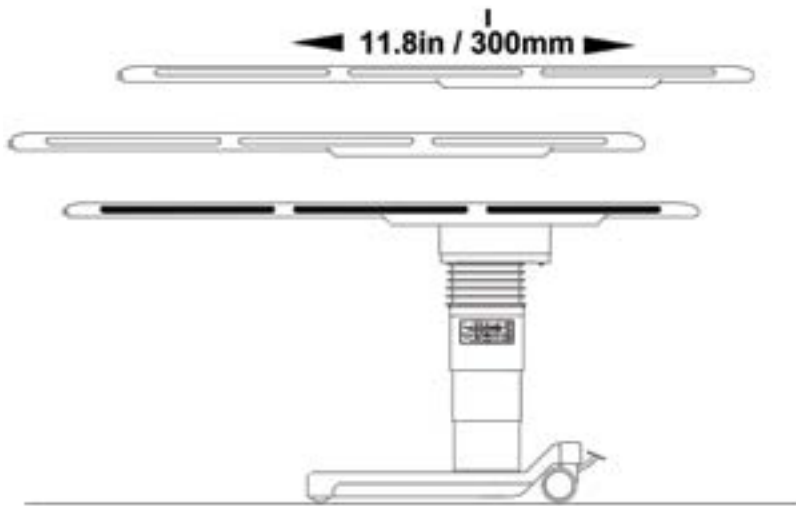