

Interventional / Endovascular Imaging

ISR

Table and Accessories Catalog



idi  
imagediagnostics.com

INTELLIGENT DESIGN ► TANGIBLE OUTCOMES



## ISR Table Description:

Mobile imaging table with motorized elevation, Trendelenburg tilt, iso-roll and manual 4-way tabletop float. For advanced vascular surgery and interventional procedures.

## ISR Table #A100-2603 (U.S. accessory rails )

### #A100-2834 (Metric accessory rails)

Includes short, recessed side rails 5" (127mm) in length at imaging end of table that facilitate use of a removable headrest extension (option) and also facilitates attachment of a catheter tray extension (option) at the narrow end of the table. Two different catheter tray extensions are available: One for the pedestal end of table and one for the imaging end of table. Table has battery backup power, 600 lb. patient load, steer assist casters, 3.75 (95 mm) tongue height, 32" (813 mm) longitudinal tabletop float, and 33" (838 mm) to 45" (1143 mm) tabletop elevation.

Includes the following:

- A100-3734: Wireless pendant hand control
- X100-2310: Tabletop Pad, 2" (51 mm) thick
- A270-0206: Wireless panning hand control
- A100-3731: Patient restraint straps (2 straps/1 pad)



Tabletop 22" (559mm) x 24" (610mm) x 84" (2134mm)

## ISR Table preferred accessory packages

### ISR Practice Starter Kit

(Peripherals - head / imaging end tray) # A100-2700

Optional package includes:

1. A100-2322: Catheter Tray Tabletop Extension – **Mounts to head / imaging end of table**
2. A100-1007: Clamp-on 8.5" US Accessory Rail
3. A310-056: Armboard, quick release, rail mounted, includes pad and strap
4. A100-2244: VAB Carbon Fiber Arm Board w / pad (Fistula Applications)
5. C000-1181: Medical Grade Anti-Fatigue Mat for maximum comfort & support
6. A610-051: Radiation Shield, Articulating, rail mounted
7. A100-3285: IV Pole, Accessory Rail Mounted, Stainless Steel, 4 bottle hooks and rail clamp

### ISR Practice Starter Kit - Cardiology # A100-2975

Optional package includes:

1. A100-1769: Catheter Tray Tabletop Extension - **Mounts to foot /pedestal end of table**
2. A100-1007: Clamp-on 8.5" US accessory rail
3. A310-056: Armboard: quick-release, rail mounted, includes pad and strap
4. A100-2244: VAB Carbon Fiber Arm Board with pad (Fistula Applications)
5. C000-1181: Medical Grade Anti-Fatigue Mat for maximum comfort & support
6. A610-051: Radiation Shield, articulating, rail-mounted
7. A100-3285: IV Pole, accessory rail-mounted, stainless steel, 4 bottle hooks and rail clamp

### ISR Accessory Kit - Vascular # A100-2977

Optional package includes:


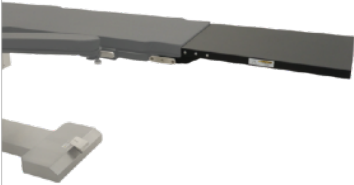



1. A100-2322: Catheter Tray Tabletop Extension – **Mounts to head / imaging end of table**
2. A100-1007: Clamp-on 8.5" US accessory rail
3. A310-056: Arm board: quick-release, rail mounted, includes pad and strap
4. A100-2244: VAB Carbon Fiber Arm Board with pad (Fistula applications)

### ISR Accessory Kit - Cardiology # A100-2978

Optional package includes:

1. A100-1769: Catheter Tray Tabletop Extension – **Mounts to foot / pedestal end of table**
2. A100-1007: Clamp-on 8.5" US Accessory Rail
3. A310-056: Armboard – rail mount
4. A100-2244: VAB Carbon Fiber Arm Board w / pad (Fistula Applications)


ISR Table Optional and Replacement Accessories:		
Photo	Description <span style="color: red;">Red text = Metric accessory rails</span>	Catalog Number
	<b>Replacement Tabletop Pad</b> 2" (51mm) thick foam pad. 19.5" (495 mm) / 22" (559 mm) / 24" (610 mm) wide x 84" (2134 mm) long.	X100-2310
	<b>Clamp-on Accessory Rail</b> Rail is 8.5" long x .375" thick x 1.125" high. <b>Clamp-on Accessory Rail</b> <span style="color: red;">10 mm x 25 mm x 216 mm clamp-on accessory rail.</span>	A100-1007 <span style="color: red;">A100-2379</span>
	<b>Armboard</b> Quick-release, rail mounted (requires Clamp-on Accessory Rail A100-1007 or <span style="color: red;">A100-2379</span> ).	A310-056 <span style="color: red;">A310-062</span>
	<b>Replacement pad for A310-056 or A310-062 Armboard</b>	C100-3329
	<b>Single use soft armband straps, package of 12</b>	C000-1555
	<b>VAB Carbon Arm Board w / Pad</b> Carbon fiber, one-piece, extra wide armboard for vascular access and fistula application procedures. Body section slides under patient/ table pad. Includes strap and pad. Clamp-on rail A100-1007 recommended.	A100-2244
	<b>Replacement Pad for A100-2244 VAB Armboard</b>	X100-2243
	<b>Armboard</b> Carbon fiber, shoulder mount (attaches by sliding thin section under patient / table pad), flat. Includes pad and strap.	A310-059
	<b>Replacement Pad for A310-059 Armboard</b>	X100-3244
	<b>Curve Spine Frame</b> Provides maximum lordosis for laminectomy, decompression, disc surgery and microdiscectomy procedures. Includes hand crank. Patient supports are adjustable to fit different patient sizes.	C000-1750
	<b>Headrest Extension</b> Articulating carbon fiber headrest. Attaches to head end of table. Includes built in anesthesia screen holder sockets. 23.125" (587 mm) wide x 12" (305 mm) long. <b>Cannot be used on rectangular tabletop.</b>	A100-2317 <span style="color: red;">A100-2505</span>

ISR Table Optional and Replacement Accessories:		
Photo	Description <span style="color: red;">Red text = Metric accessory rails</span>	Catalog Number
	<b>Replacement Pad for A100-2317 or A100-2505 Headrest Extension</b>	X100-2325
	<b>Peripheral Extension</b> Carbon Fiber 12" (305 mm) table extension for imaging end of table. Table top must have mini rails as shown.	A100-2486
	<b>Replacement Pad for A100-2486 Peripheral Extension</b>	X100-2503
	<b>Tabletop Catheter Tray Extension</b> 20.75" (527 mm) x 27" (686 mm) long. <b>Mounts to imaging end of table only.</b> The maximum weight this extension can support is 20 lb or 9 kg.	A100-2322
	<b>Radiation Shield - .5mm lead equivalence</b> Articulating, rail-mounted. Attaches to table side-rail at base end of table, either side, adjustable. Articulating joint allows shield to conform to table contour. Snap-on patient side upper section is removable. 41.375" (1051mm) wide x 32.125" (816 mm) high x 2.625" (16 mm) deep.	A610-051 <span style="color: red;">A610-106</span>
	<b>Compass Contour Radiation Shield</b> Contoured ceiling suspended acrylic barrier radiation protection with a lead equivalence of .5mm. Enhance eye, head and neck protection from x-ray scatter. Shielding area is ~30" (76 cm) x 23.5" (60 cm) with a 10" (254 mm) radius along the bottom edge of the leaded acrylic shield. Does not include ceiling suspension arm.	A610-0200
	<b>Compass Frame Radiation Shield</b> Rectangular ceiling suspended acrylic barrier radiation protection with a lead equivalence of 0.5mm. Enhance eye, head and neck protection from x-ray scatter. Shielding area is ~ 30" (76 cm) x 23.5" (60 cm). Does not include ceiling suspension arm.	A610-0215

ISR Table Optional and Replacement Accessories:		
Photo	Description <span style="color: red;">Red text = Metric accessory rails</span>	Catalog Number
	<p><b>Compass Contour and Frame Drape</b>            Custom shaped drape for the Compass Contour and Frame Radiation Shield. Slips over the entire shield and is held in place with an elastic band. Sterile. Case of 20.</p>	C610-0256
	<p><b>Protego™ Radiation Shield System</b>            Mobile Imaging Shielding System for advanced cardiovascular procedures. The system when properly installed and monitored may provide full body X-ray protection for all members of the clinical team. The system mounts to existing table rails and is equipped with a mobile support stand supporting the top shield. There is a floor shield providing additional side coverage, in addition to patient sterile and non-sterile drapes and shields.</p>	Please consult with Image Diagnostics to configure correct Protego system options
	<p><b>Anesthesia Screen Holder</b>            Flexible "hoop" style. Mounts directly-A100-2317 Headrest Extension. For use on accessory rails - a pair of rotating side rail clamps, part number C000-0746, are required. Rotating side rail clamps and clamp-on accessory rails are not required if table is ordered with the A100-2317 Headrest Extension.</p>	C000-0965
	<p><b>ISR Wireless Panning Hand Control</b>            Included with table</p>	A270-0206
	<p><b>Wireless ISR Pendant Hand Control</b>            Included with table</p>	A100-3734

ISR Table Optional and Replacement Accessories:		
Photo	Description <span style="color: red;">Red text = Metric accessory rails</span>	Catalog Number
	<b>Heavy Duty Medical Grade Anti-Fatigue Mat</b> 24" (610 mm) x 66" (1676 mm) x 3/4" (19 mm) IDI heavy-duty anti-fatigue mat. One-piece construction for high-traffic areas. Safe & non-toxic, contains no plasticizers and will not leach over time. Durable brushed texture.	C000-1181
	<b>Simple clamp</b> Accommodates a round accessory post. All metal design clamp mounts and locks anywhere on standard OR table rail. Autoclavable.	C100-3046
	<b>Snap-on Blade Clamp</b> Accommodates a blade-style accessory post. Quickly snaps onto rail, push button to release clamp. Mounts and locks anywhere on standard OR table rail. Prevents false rail engagement. Autoclavable.	C100-3043
	<b>Simple Blade Clamp</b> Accommodates blade-style accessory post. Mounts and locks anywhere on standard OR table rail ends or at rail notches. A cost-effective solution.	C000-1115
	<b>Snap-on Socket</b> Accommodates a round accessory post. Quickly snaps onto rail and pull bottom lever to release. Mounts and locks anywhere on standard OR table rail. Lithotomy wheel for accessory angle adjustments.	C000-0746

### ISR Table Optional and Replacement Accessories:

Photo	Description <span style="color: red;">Red text = Metric accessory rails</span>	Catalog Number
	<b>Restraint Straps</b> Set of two patient restraint straps with buckle locking mechanism and one pad for patient safety and comfort	A100-3731

### Shipping Crate Specifications

- Packaged weight (approx.): 890 lbs. (404 kg)
- Length: 96" (2438 mm).
- Width: 42" (1067 mm).
- Height: 47" (1194 mm)





Since 1986, IDI has been a manufacturer of innovative products for surgical, interventional and diagnostic imaging.

Our product focus is in electro-mechanical solutions for fluoroscopic and digital radiographic applications.

IDI is committed to superior product design, value and customer support in everything we do.

**Image Diagnostics, Inc.**  
310 Authority Drive  
Fitchburg, Massachusetts  
01420

For customer service:  
[sales@imagediagnostics.com](mailto:sales@imagediagnostics.com)  
Telephone: 978.829.0009

**INTELLIGENT DESIGN ► TANGIBLE OUTCOMES**



ORDER: [sales@imagediagnostics.com](mailto:sales@imagediagnostics.com)  
SUPPORT: [techsupport@imagediagnostics.com](mailto:techsupport@imagediagnostics.com)  
WEBSITE: <https://imagediagnostics.com>

+1-978-829-0009  
Monday through Friday

