


100UC+
UROLOGY TABLE

100US
GYN / FERTILITY TABLE

Tables and Accessories Catalog

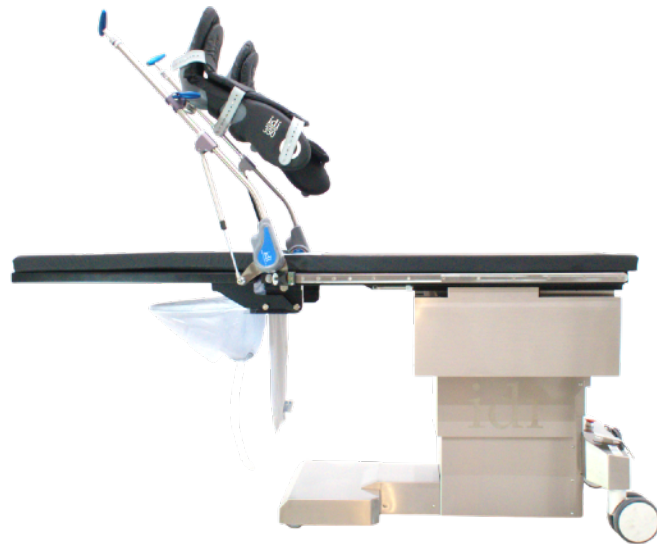

100UC+ Urology Table
shown with optional ComfortFin 350 leg holders



100US GYN / Fertility Table


INTELLIGENT DESIGN ► **TANGIBLE OUTCOMES**

idi
imagediagnostics.com

100UC+ Urology/Multi-purpose Table



100UC+ shown with optional ComfortFin® stirrups (#C100-3042)

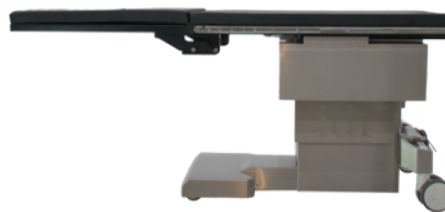

100UC+ shown with lightweight aluminum boarding extension (Included)

100UC+ description:

Imaging table for urology and multi-purpose use. Urological accessories are easily removable so the table can be quickly configured for multi-purpose C-arm imaging or general patient positioning for non-urological procedures.

100UC+ #A100-3099 (120VAC with U.S. Accessory Rails)

Includes the following:

- X100-1818: Tabletop pad, 2" (51mm) thick
- A100-3543: Aluminum boarding extension
- X100-1817: Tabletop extension pad
- A100-3733: Wireless Pendant hand control
- A100-3715: Wireless Footswitch control
- U100-3025: Removable x-ray shield flap, for perineal end of table
- A100-3883: Patient restraint strap with pad
- A310-056: Arm board, rail mounted, with pad and strap
- Z100-3359: Urology drain bag support hoop
- C000-0593: Disposable drain bag & hose assembly, sterile, box of 20.
- C000-0492: Disposable covers for foot switch, box of 50
- C000-0645: Disposable tabletop pad protector (qty of 1)
- C000-0612: 5 gallon collection bag (qty of 1)

100UC+ Package Special #A100-4114

Includes the following:

- A100-3099: Aspect 100UC+ Urology Table with standard accessories
- A100-2246: Two ComfortFin® 350 stirrups with clamps for US rails
- A310-056: One rail-mounted arm board
- C000-0593: Disposable drain bag & hose assembly, sterile, box of 20
- C100-3044: One ComfortFin stirrup storage cart

100UC+ Package Special #A100-4115

Includes the following:

- A100-3099: Aspect 100UC+ Urology Table with standard accessories
- A100-2247: Two ComfortFin 500 stirrups with clamps for US rails
- A310-056: One rail-mounted arm board
- C000-0593: Disposable drain bag & hose assembly, sterile, box of 20
- C100-3044: One ComfortFin stirrup storage cart

Photo	Description	Catalog Number
	IDI ComfortFin 350 Lithotomy Leg holders (Stirrups pair only) ComfortFin Lithotomy Leg holders are designed to improve infection control with durable, easy to clean silicone straps. Includes integrated lateral fin to help protect the head of the fibula and the peroneal nerve. 4" more working space with an increased range of motion. Patient weight capacity: 350 lbs (159 kg).	C100-3042
	Replacement Boot with Straps for IDI C100-3042 ComfortFin 350 Lithotomy Leg Holders	C100-3080
	Kit: IDI ComfortFin350 Lithotomy Leg holders C100-3042 with C100-3046 economy side-rail clamps	A100-2246
	IDI ComfortFin 500 Lithotomy Leg holders (Stirrups pair only) ComfortFin Lithotomy Leg holders are designed to improve infection control with durable, easy to clean silicone straps. Includes integrated lateral fin to help protect the head of the fibula and the peroneal nerve. Boots are sized for today's heavier patients. 4" (102mm) more working space with an increased range of motion. Patient weight capacity: 500 lbs. (227 kg).	C100-3045
	Replacement Boot with Straps for IDI C100-3045 ComfortFin 500 Lithotomy Leg Holders	C100-3079
	Kit: IDI ComfortFin 500 Lithotomy Leg holders C100-3045 with C100-3046 economy side-rail clamps.	A100-2247
	Lithotomy Leg Holder Cart Easy to push/pull with ergonomic handles. Stirrups are isolated to prevent damage. Clamp storage compartment. Designed to safely store ComfortFin stirrups when not in use.	C100-3044
	Replacement Tabletop Pad 2" (51mm) thick foam pad. 28" (711mm) wide x 46 1/2" (1181mm) long. Replacement Extension Tabletop Pad 2 " (51mm) thick foam pad. 28" (71mm) wide x 34" (864mm) long.	C100-3070
	Single use soft armband straps, package of 12	C000-1555

Photo	Description	Catalog Number
	Split Leg Holders For use with Aspect 100UC + Tables. Includes pads. Prone percutaneous, lateral or modified lateral nephrostomy. Includes two split leg holder clamps # C000-1115	C000-1146
	Right replacement pad for Spilt Leg Holder C000-1146	C000-1357
	Left replacement pad for Spilt Leg Holder C000-1146	C000-1358
	Replacement Split Leg Holder Clamp For use with split leg holders. # C000-1146	C000-1115
	Simple clamp Accommodates a round accessory post. All metal clamp mounts and locks anywhere on a standard OR table rail. Autoclaveable.	C100-3046
	Snap-on Blade Clamp Accommodates a blade-style accessory post. Quickly snaps onto rail, push button to release clamp. Mounts and locks anywhere on a standard OR table rail. Prevents false rail engagement. Autoclaveable.	C100-3043
	Simple Blade Clamp Accommodates blade-style accessory post. Mounts and locks anywhere on standard OR table rail ends or at rail notches. A cost-effective solution.	C000-1115
	Snap-on Socket Accommodates a round accessory post. Quickly snaps onto rail and pull bottom lever to release. Mounts and locks anywhere on standard OR table rail. Lithotomy wheel for accessory angle adjustments.	C000-0746



Photo	Description	Catalog Number
	Armboard Quick-release, rail mounted (requires Clamp-on Accessory Rail A100-1007 or A100-2379 for use on tabletop extension.	A310-056
	Replacement pad for A310-056 or A310-062 Armboard	C100-3329
	Clamp-on Accessory Rail (For use on tabletop extension) Rail is 8 1/2" long x 3/8" thick x 1-1/8" high.	A100-2218
	Stirrups Mounting Rail extension Kit Included with table. For use with leg-holders (Leg-holders are sold separately).	A100-1876
	Wireless Hand Control for 100UC+ Controls height up/down, trendelenburg/reverse, slide, cross table movements. For tables manufactured after July 2025. Included with table.	A100-3733
	Wireless Foot Switch Control for 100UC+ Included with table. For tables manufactured after July 2025.	A100-3715
	Hand Control for 100UCPlus and 100UC Controls height up/down, trendelenburg/reverse, slide, cross table movements. For tables manufactured before July 2025 Included with table.	K000-0158
	Foot Switch Control for 100UCPlus. For tables manufactured before July 2025 Included with table.	A100-1540





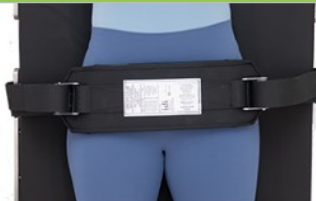


Photo	Description	Catalog Number
	Compass Frame Radiation Shield Rectangular ceiling suspended acrylic barrier radiation protection with a lead equivalence of 0.5mm. Enhance eye, head and neck protection from x-ray scatter. Shielding area is ~30" (760 mm) x 23.5" (600 mm). Does not include ceiling suspension arm.	A610-0215
	Compass Contour and Frame Drape Custom shaped drape for the Compass Contour and Frame Radiation Shield. Slips over the entire shield and is held in place with an elastic band. Sterile. Case of 20.	C610-0256
	IV Pole Stainless steel with four bottle hooks Range: 42" (1070 mm) / 65" (1660 mm) Tube diameter: .98" (25 mm) / 0.079" (18 mm) IV Pole Clamp C000-1325 is included.	A100-3285
	Vinyl Urology Table Covers Box of 20.	C000-0645
	Uro Drain Bag with Hose Disposable drain bag and hose assembly. Box of 20. For tables manufactured after July 2025.	Sterile: C000-0593 Non-sterile: C000-1111
	Disposable Covers for Aspect 100UC+ Foot Controls 24" (610 mm) x 20" (508 mm) Box of 20.	C000-0492
	Uro Drain Support Hoop Collapsible & removable drain bag support hoop. For tables manufactured before July 2025	C000-0595

Photo	Description	Catalog Number
	Restraint Strap Patient restraint strap with hooks to attach to accessory rails on each side of the tabletop, a pad for patient safety and comfort and buckles to adjust to the patient girth	A100-3883
	Uro Collection Unit 5 gallon/18.9 liter capacity, Box of 10. Soft sided Vinyl container gently expands into a cube as it collects waste fluid when floor drain or suction is unavailable. Internal baffle system enables container to fill evenly as powder converts liquid to gel. Built in zipper and handles allow for storage, closure of drain bag and hose assembly for enhanced infection control during transport and disposal.	C000-0612
	Uro Collection Unit 15 liter capacity, Box of 10. Flexible plastic sided container expands into a cube as it collects waste fluid when floor drain or suction is unavailable. Internal baffle system enables container to fill evenly as powder converts liquid to gel. Integrated valve system allows user to shut off the flow of liquid from the drain bag in order to change collection units.	C000-1307
100 UC+ Shipping Crate Specifications		
<ul style="list-style-type: none"> • Weight (approx.): 725 lbs. (330 kg) • Length: 96" (2438 mm) • Width: 42" (1067 mm) • Height: 47" (1194 mm) 		

100US Gyn/Fertility Clinic Table




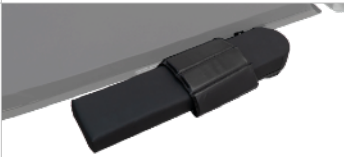

100US A100-1479

GYN/Fertility Mobile Imaging Table

An economical C-arm table for hysterosalpingograms

Includes the following:

- A100-758: Foot control
- X100-1687: Tabletop pad, 24" x 46" (610 mm x 1168 mm) 1" (25.4 mm) thick
- A100-3883: Patient restraint strap
- A100-2380: Leg support tabletop extension with pad
- A100-2248: Basic Knee Crutches [pair] with a pair of rotating side rail clamps [C000-0746]

	Replacement Tabletop Pad 2" (51 mm) thick foam pad. 24" (610 mm) wide x 46" (1168 mm) long.	X100-1687
	Armboard Quick-release, rail mounted.	A310-056
	Replacement pad for A310-056 or A310-062 Armboard	C100-3288
	Single use soft armband straps, package of 12	C000-1555
	Hand control for 100US	K000-0041

	<p>Basic Knee Crutches For patients up to 400 lbs. (180 kg) Basic pair with pads. Includes two C000-0746 rotating side rail clamps.</p>	<p>A100-2248</p>
	<p>Foot Control 100US (Included with table.)</p>	<p>A100-758</p>
	<p>IV Pole Stainless steel with four bottle hooks Range: 42" (1070 mm) / 65" (1660) mm Tube diameter: .98" (25 mm) / .079" (18 mm) IV Pole Clamp C000-1325 is included.</p>	<p>C000-1324</p>
	<p>Compass Contour Radiation Shield Contoured ceiling suspended acrylic barrier radiation protection with a lead equivalence of .5mm. Enhance eye, head and neck protection from x-ray scatter. Shielding area is ~30" (760 mm) x 23.5" (600 mm) with a 10" (254 mm) radius along the bottom edge of the leaded acrylic shield. Does not include ceiling suspension arm.</p>	<p>A610-0200</p>
	<p>Compass Frame Radiation Shield Rectangular ceiling suspended acrylic barrier radiation protection with a lead equivalence of 0.5mm. Enhance eye, head and neck protection from x-ray scatter. Shielding area is ~30" (760 mm) x 23.5" (600 mm). Does not include ceiling suspension arm.</p>	<p>A610-0215</p>
	<p>Compass Contour and Frame Drape Custom shaped drape for the Compass Contour and Frame Radiation Shield. Slips over the entire shield and is held in place with an elastic band. Sterile. Case of 20.</p>	<p>C610-0256</p>
	<p>Restraint Strap Patient restraint strap with hooks to attach to accessory rails on each side of the tabletop, a pad for patient safety and comfort and buckles to adjust to the patient girth</p>	<p>A100-3883</p>

	Simple clamp Accommodates a round accessory post. All metal design clamp mounts and locks anywhere on standard OR table rail. Autoclavable.	C100-3046
	Snap-on Blade Clamp Accommodates a blade-style accessory post. Quickly snaps onto rail, push button to release clamp. Mounts and locks anywhere on standard OR table rail. Prevents false rail engagement. Autoclavable.	C100-3043
	Simple Blade Clamp Accommodates blade-style accessory post. Mounts and locks anywhere on standard OR table rail ends or at rail notches. A cost-effective solution.	C000-1115
	Snap-on Socket Accommodates a round accessory post. Quickly snaps onto rail and pull bottom lever to release. Mounts and locks anywhere on standard OR table rail. Lithotomy wheel for accessory angle adjustments.	C000-0746
	Uro Collection Unit 5 gallon/18.9 liter capacity capacity, Box of 10. Soft sided Vinyl container gently expands into a cube as it collects waste fluid when floor drain or suction is unavailable. Internal baffle system enables container to fill evenly as powder converts liquid to gel. Built in zipper and handles allow for storage, closure of drain bag and hose assembly for enhanced infection control during transport and disposal.	C000-0612

100US Shipping Crate Specifications

- Weight (approx.): 815 lbs. (370 kg)
- Length: 96" (2438 mm)
- Width: 42" (1067 mm)
- Height: 47" (1194 mm)

Aspect Urology & GYN Table Specifications & Accessories

	100 UC+	100 US
Elevate (Height with pad)	29.25" – 41" (74.3 cm – 104.14 cm)	27" – 39" (686 mm – 991 mm)
Trendelenburg Tilt	± 12°	NA
Lateral travel	± 3.5", 7" total (± 8.8 cm, 17.6cm total)	NA
Longitudinal Tabletop Motion	13.75" (35cm)	NA
Tabletop Size	28"W x 46.5"L (71cm 118cm)	24"W x 46"L (610 mm x 1168 mm)
Radiolucent Tabletop Extension Size	28/24"W x 34"L (71/61cm x 86cm)	NA
Leg Boarding Tabletop Extension	28"W x 32"L (71cm 81cm)	
Imaging Area	Main tabletop: 29"W x 23"L (73.5cmmm x 58cm)	22"W x 32"L (559 mm x 813 mm)
Imaging Area with optional radiolucent table extension	54" (137 cm)	NA
Patient Weight Capacity	Upper body on main table: 550lb. (249 kg) without tabletop extensions and accessories	450 lbs (204 kg)
Table Weight	~724 lb. (~329 kg)	375 lbs ((170 kg)
Multi-Caster/3 Postion Multilock system	Total Lock, Steer & Incline	NA
Tabletop Pad 2" (51mm) thick)	•	•
Tabletop extension with pad	•	For leg support only
Hand control pendant	•	•
Footswitch control	•	◇
Patient restraint straps	•	•
Removeable drain bag support hoop	•	NA
Disposable drain bag & hose assembly, box of 20	•	NA
Removable radiation shield flap for perineal end of table	•	NA
Armboard, rail-mounted with pad and strap	•	◇
Split Leg Stirrups	◇	NA
Disposable covers for footswitch, box of 50	•	◇
Bluetooth compatible foot and hand operator controls	•	NA
Emergency backup battery	•	•
Emergency stop pushbutton	•	•
	•=Included ◇=Optional NA=Not available	



Since 1986, IDI has been a manufacturer of innovative products for surgical, interventional and diagnostic imaging.

Our product focus is in electro-mechanical solutions for fluoroscopic and digital radiographic applications.

IDI is committed to superior product design, value and customer support in everything we do.

Image Diagnostics, Inc.
310 Authority Drive
Fitchburg, Massachusetts
01420

For customer service:
sales@imagediagnostics.com
Telephone: 978.829.0009

INTELLIGENT DESIGN ► TANGIBLE OUTCOMES



ORDER: sales@imagediagnostics.com
SUPPORT: techsupport@imagediagnostics.com
WEBSITE: <https://imagediagnostics.com>

+1-978-829-0009
Monday through Friday



This page left intentionally blank